

Factor IX19. GlycoPEGylation of Factor IX produced in CHO cells

This example sets forth the preparation of asialoFactor IX and its sialylation with CMP-sialic acid-PEG.

- 5       **Desialylation of rFactor IX.** A recombinant form of Coagulation Factor IX (rFactor IX) was made in CHO cells. 6000 IU of rFactor IX were dissolved in a total of 12 mL USP H<sub>2</sub>O. This solution was transferred to a Centricon Plus 20, PL-10 centrifugal filter with another 6 mL USP H<sub>2</sub>O. The solution was concentrated to 2 mL and then diluted with 15 mL 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl<sub>2</sub>, 0.05% NaN<sub>3</sub> and then reconcentrated.
- 10      The dilution/concentration was repeated 4 times to effectively change the buffer to a final volume of 3.0 mL. Of this solution, 2.9 mL (about 29 mg of rFactor IX) was transferred to a small plastic tube and to it was added 530 mU  $\alpha$ 2-3,6,8-Neuraminidase-agarose conjugate (*Vibrio cholerae*, Calbiochem, 450  $\mu$ L). The reaction mixture was rotated gently for 26.5 hours at 32 °C. The mixture was centrifuged 2 minutes at 10,000 rpm and the supernatant
- 15      was collected. The agarose beads (containing neuraminidase) were washed 6 times with 0.5 mL 50 mM Tris-HCl pH 7.12, 1 M NaCl, 0.05% NaN<sub>3</sub>. The pooled washings and supernatants were centrifuged again for 2 minutes at 10,000 rpm to remove any residual agarose resin. The pooled, desialylated protein solution was diluted to 19 mL with the same buffer and concentrated down to ~2 mL in a Centricon Plus 20 PL-10 centrifugal filter. The
- 20      solution was twice diluted with 15 mL of 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 0.05% NaN<sub>3</sub> and reconcentrated to 2 mL. The final desialylated rFactor IX solution was diluted to 3 mL final volume (~10 mg/mL) with the Tris Buffer. Native and desialylated rFactor IX samples were analyzed by IEF-Electrophoresis. Isoelectric Focusing Gels (pH 3-7) were run using 1.5  $\mu$ L (15  $\mu$ g) samples first diluted with 10  $\mu$ L Tris buffer and mixed with 12  $\mu$ L
- 25      sample loading buffer. Gels were loaded, run and fixed using standard procedures. Gels were stained with Colloidal Blue Stain (Figure 154), showing a band for desialylated Factor IX.

- Preparation of PEG (1 kDa and 10 kDa)-SA-Factor IX.** Desialylated rFactor-IX (29 mg, 3 mL) was divided into two 1.5 mL (14.5 mg) samples in two 15 mL centrifuge
- 30      tubes. Each solution was diluted with 12.67 mL 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 0.05% NaN<sub>3</sub> and either CMP-SA-PEG-1k or 10k (7.25  $\mu$ mol) was added. The tubes were

inverted gently to mix and 2.9 U ST3Gal3 (326  $\mu$ L) was added (total volume 14.5 mL). The tubes were inverted again and rotated gently for 65 hours at 32 °C. The reactions were stopped by freezing at -20 °C. 10  $\mu$ g samples of the reactions were analyzed by SDS-PAGE. The PEGylated proteins were purified on a Toso Haas Biosep G3000SW (21.5 x 30 cm, 13  $\mu$ m) HPLC column with Dulbecco's Phosphate Buffered Saline, pH 7.1 (Gibco), 6 mL/min. The reaction and purification were monitored using SDS Page and IEF gels. Novex Tris-Glycine 4-20% 1 mm gels were loaded with 10  $\mu$ L (10  $\mu$ g) of samples after dilution with 2  $\mu$ L of 50 mM Tris-HCl, pH 7.4, 150 mM NaCl, 0.05% NaN<sub>3</sub> buffer and mixing with 12  $\mu$ L sample loading buffer and 1  $\mu$ L 0.5 M DTT and heated for 6 minutes at 85 °C. Gels were stained with Colloidal Blue Stain (Figure 155) showing a band for PEG (1 kDa and 10 kDa)-SA-Factor IX.

#### 20. Direct Sialyl-GlycoPEGylation of Factor IX

This example sets forth the preparation of sialyl-PEGylation of Factor IX without prior sialidase treatment.

**Sialyl-PEGylation of Factor-IX with CMP-SA-PEG-(10 kDa).** Factor IX (1100 IU), which was expressed in CHO cells and was fully sialylated, was dissolved in 5 mL of 20 mM histidine, 520 mM glycine, 2% sucrose, 0.05% NaN<sub>3</sub> and 0.01% polysorbate 80, pH 5.0. The CMP-SA-PEG-(10 kDa) (27 mg, 2.5  $\mu$ mol) was then dissolved in the solution and 1 U of ST3Gal3 was added. The reaction was complete after gently mixing for 28 hours at 32°C. The reaction was analyzed by SDS-PAGE as described by Invitrogen. The product protein was purified on an Amersham Superdex 200 (10 x 300 mm, 13  $\mu$ m) HPLC column with phosphate buffered saline, pH 7.0 (PBS), 1 mL/min. R<sub>t</sub> = 9.5 min.

**Sialyl-PEGylation of Factor-IX with CMP-SA-PEG-(20 kDa).** Factor IX (1100 IU), which was expressed in CHO cells and was fully sialylated, was dissolved in 5 mL of 20 mM histidine, 520 mM glycine, 2% sucrose, 0.05% NaN<sub>3</sub> and 0.01% polysorbate 80, pH 5.0. The CMP-SA-PEG-(20 kDa) (50 mg, 2.3  $\mu$ mol) was then dissolved in the solution and CST-II was added. The reaction mixture was complete after gently mixing for 42 hours at 32°C. The reaction was analyzed by SDS-PAGE as described by Invitrogen.

The product protein was purified on an Amersham Superdex 200 (10 x 300 mm, 13  $\mu$ m) HPLC column with phosphate buffered saline, pH 7.0 (Fisher), 1 mL/min.  $R_t$  = 8.6 min.

#### 21. Sialic Acid Capping of GlycoPEGylated Factor IX

This examples sets forth the procedure for sialic acid capping of sialyl-glycoPEGylated peptides. Here, Factor-IX is the exemplary peptide.

**Sialic acid capping of N-linked and O-linked Glycans of Factor-IX-SA-PEG (10 kDa).** Purified r-Factor-IX-PEG (10 kDa) (2.4 mg) was concentrated in a Centricon® Plus 20 PL-10 (Millipore Corp., Bedford, MA) centrifugal filter and the buffer was changed to 50 mM Tris-HCl pH 7.2, 0.15 M NaCl, 0.05%  $\text{NaN}_3$  to a final volume of 1.85 mL. The protein solution was diluted with 372  $\mu$ L of the same Tris buffer and 7.4 mg CMP-SA (12  $\mu$ mol) was added as a solid. The solution was inverted gently to mix and 0.1 U ST3Gal1 and 0.1 U ST3Gal3 were added. The reaction mixture was rotated gently for 42 hours at 32 °C.

A 10  $\mu$ g sample of the reaction was analyzed by SDS-PAGE. Novex Tris-Glycine 4-12% 1 mm gels were performed and stained using Colloidal Blue as described by Invitrogen. Briefly, samples, 10  $\mu$ L (10  $\mu$ g), were mixed with 12  $\mu$ L sample loading buffer and 1  $\mu$ L 0.5 M DTT and heated for 6 minutes at 85 °C (Figure 156, lane 4).

#### Factor VIIa

##### 22. GlycoPEGylation of Recombinant Factor VIIa produced in BHK cells

This example sets forth the PEGylation of recombinant Factor VIIa made in BHK cells.

**Preparation of Asialo-Factor VIIa.** Recombinant Factor VIIa was produced in BHK cells (baby hamster kidney cells). Factor VIIa (14.2 mg) was dissolved at 1 mg/ml in buffer solution (pH 7.4, 0.05 M Tris, 0.15 M NaCl, 0.001 M  $\text{CaCl}_2$ , 0.05%  $\text{NaN}_3$ ) and was incubated with 300 mU/mL sialidase (*Vibrio cholera*)-agarose conjugate for 3 days at 32 °C. To monitor the reaction a small aliquot of the reaction was diluted with the appropriate buffer and an IEF gel performed according to Invitrogen procedures (Figure 157). The mixture was centrifuged at 3,500 rpm and the supernatant was collected. The resin was washed three times (3 $\times$ 2 mL) with the above buffer solution ( pH 7.4, 0.05 M Tris, 0.15 M NaCl, 0.05%  $\text{NaN}_3$ ) and the combined washes were concentrated in a Centricon-Plus-20. The remaining

solution was buffer exchanged with 0.05 M Tris (pH 7.4), 0.15 M NaCl, 0.05% NaN<sub>3</sub> to a final volume of 14.4 mL.

**Preparation of Factor VIIa-SA-PEG (1 kDa and 10 kDa).** The desialylation of Factor VIIa solution was split into two equal 7.2 ml samples. To each sample was added either CMP-SA-5-PEG(1 kDa) (7.4 mg) or CMP-SA-5-PEG(10 kDa) (7.4 mg). ST3Gal3 (1.58U) was added to both tubes and the reaction mixtures were incubated at 32°C for 96 hrs. The reaction was monitored by SDS-PAGE gel using reagents and conditions described by Invitrogen. When the reaction was complete, the reaction mixture was purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The combined fractions containing the product were concentrated at 4°C in Centricon-Plus-20 centrifugal filters (Millipore, Bedford, MA) and the concentrated solution reformulated to yield 1.97 mg (bicinchoninic acid protein assay, BCA assay, Sigma-Aldrich, St. Louis MO) of Factor VIIa-PEG. The product of the reaction was analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples were dialyzed against water and analyzed by MALDI-TOF. Figure 158 shows the MALDI results for native Factor VIIa. Figure 159 contains the MALDI results for Factor VIIa PEGylated with 1 kDa PEG where peak of Factor VIIa PEGylated with 1kDa PEG is evident. Figure 160 contains the MALDI results for Factor VIIa PEGylated with 10 kDa PEG where a peak for Factor VIIa PEGylated with 10 kDa PEG is evident. Figure 161 depicts the SDS-PAGE analysis of all of the reaction products, where a band for Factor VIIa-SA-PEG (10 kDa) is evident.

### Follicle Stimulating Hormone (FSH)

#### 23. GlycoPEGylation of human pituitary-derived FSH

This example illustrates the assembly of a conjugate of the invention. Follicle Stimulating Hormone (FSH) is desialylated and then conjugated with CMP-(sialic acid)-PEG.

**Desialylation of Follicle Stimulating Hormone.** Follicle Stimulating Hormone (FSH) (Human Pituitary, Calbiochem Cat No. 869001), 1 mg, was dissolved in 500 µL 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl<sub>2</sub>. This solution, 375 µL, was transferred to a small plastic tube and to it was added 263 mU Neuraminidase II (*Vibrio cholerae*). The reaction mixture was shaken gently for 15 hours at 32 °C. The reaction mixture was added to

N-(*p*-aminophenyl)oxamic acid-agarose conjugate, 600  $\mu$ L, pre-equilibrated with 50 mM Tris-HCl pH 7.4, 150 mM NaCl and 0.05% NaN<sub>3</sub> and gently rotated 6.5 hours at 4 °C. The suspension was centrifuged for 2 minutes at 14,000 rpm and the supernatant was collected. The beads were washed 5 times with 0.5 mL of the buffer and all supernatants were pooled.

The enzyme solution was dialyzed (7000 MWCO) for 15 hours at 4 °C with 2 L of a solution containing 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub>, and then twice for 4 hours at 4 °C into 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub>. The solution was concentrated to 2  $\mu$ g/ $\mu$ L by Speed Vac and stored at -20 °C. Reaction samples were analyzed by IEF gels (pH 3-7) (Invitrogen) (Figure 162).

#### Preparation of human pituitary-derived SA-FSH and PEG-SA-Follicle

**Stimulating Hormone.** Desialylated FSH (100  $\mu$ g, 50  $\mu$ L) and CMP-sialic acid or CMP-SA-PEG (1 kDa or 10 kDa) (0.05  $\mu$ mol) were dissolved in 13.5  $\mu$ L H<sub>2</sub>O (adjusted to pH 8 with NaOH) in 0.5 mL plastic tubes. The tubes were vortexed briefly and 40 mU ST3Gal3 (36.5  $\mu$ L) was added (total volume 100  $\mu$ L). The tubes were vortexed again and shaken gently for 24 hours at 32 °C. The reactions were stopped by freezing at -80 °C. Reaction samples of 15  $\mu$ g were analyzed by SDS-PAGE (Figure 163), IEF gels (Figure 164) and MALDI-TOF. Native FSH was also analyzed by SDS-PAGE (Figure 165)

**Analysis of SDS PAGE and IEF Gels of Reaction Products.** Novex Tris-Glycine 8-16% 1 mm gels for SDS PAGE analysis were purchased from Invitrogen. 7.5  $\mu$ L (15  $\mu$ g) of FSH reaction samples were diluted with 5  $\mu$ L of 50 mM Tris-HCl, pH 7.4, 150 mM NaCl, 0.05% NaN<sub>3</sub> buffer, mixed with 15  $\mu$ L sample loading buffer and 1  $\mu$ L 9 M  $\mu$ -mercaptoethanol and heated for 6 minutes at 85 °C. Gels were run as directed by Invitrogen and stained with Colloidal Blue Stain (Invitrogen).

FSH samples (15  $\mu$ g) were diluted with 5  $\mu$ L Tris buffer and mixed with 15  $\mu$ L sample loading buffer (Figure 162). The samples were then applied to Isoelectric Focusing Gels (pH 3-7) (Invitrogen) (Figure 165). Gels were run and fixed as directed by Invitrogen and then stained with Colloidal Blue Stain.

24. GlycoPEGylation of recombinant FSH produced recombinantly in CHO cells

This example illustrates the assembly of a conjugate of the invention. Desialylated FSH was conjugated with CMP-(sialic acid)-PEG.

- 5       **Preparation of recombinant Asialo-Follicle Stimulation Hormone.** Recombinant Follicle Stimulation Hormone (rFSH) produced from CHO was used in these studies. The 7,500 IU of rFSH was dissolved in 8 mL of water. The FSH solution was dialyzed in 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl<sub>2</sub> and concentrated to 500  $\mu$ L in a Centricon Plus 20 centrifugal filter. A portion of this solution (400  $\mu$ L) (~0.8 mg FSH) was transferred to a
- 10   small plastic tube and to it was added 275 mU Neuraminidase II (*Vibrio cholerae*). The reaction mixture was mixed for 16 hours at 32 °C. The reaction mixture was added to prewashed N-(p-aminophenyl)oxamic acid-agarose conjugate (800  $\mu$ L) and gently rotated for 24 hours at 4 °C. The mixture was centrifuged at 10,000 rpm and the supernatant was collected. The beads were washed 3 times with 0.6 mL Tris-EDTA buffer, once with 0.4 mL
- 15   Tris-EDTA buffer and once with 0.2 mL of the Tris-EDTA buffer and all supernatants were pooled. The supernatant was dialyzed at 4 °C against 2 L of 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub> and then twice more against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub>. The dialyzed solution was then concentrated to 420  $\mu$ L in a Centricon Plus 20 centrifugal filter and stored at -20 °C.
- 20       Native and desialylated rFSH samples were analyzed by SDS-PAGE and IEF (Figure 166). Novex Tris-Glycine 8-16% 1 mm gels were purchased from Invitrogen. Samples (7.5  $\mu$ L, 15  $\mu$ g) samples were diluted with 5  $\mu$ L of 50 mM Tris-HCl, pH 7.4, 150 mM NaCl, 0.05% NaN<sub>3</sub> buffer, mixed with 15  $\mu$ L sample loading buffer and 1  $\mu$ L 9 M  $\beta$ -mercaptoethanol and heated for 6 minutes at 85 °C. Gels were run as directed by Invitrogen
- 25   and stained with Colloidal Blue Stain (Invitrogen). Isoelectric Focusing Gels (pH 3-7) were purchased from Invitrogen. Samples (7.5  $\mu$ L, 15  $\mu$ g) were diluted with 5  $\mu$ L Tris buffer and mixed with 15  $\mu$ L sample loading buffer. Gels were loaded, run and fixed as directed by Invitrogen. Gels were stained with Colloidal Blue Stain. Samples of native and desialylated FSH were also dialyzed against water and analyzed by MALDI-TOF.
- 30       **Sialyl-PEGylation of recombinant Follicle Stimulation Hormone.** Desialylated FSH (100  $\mu$ g, 54  $\mu$ L) and CMP-SA-PEG (1 kDa or 10 kDa) (0.05  $\mu$ mol) were dissolved in 28

μL 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2 in 0.5 mL plastic tubes. The tubes were vortexed briefly and 20 mU of ST3Gal3 was added (total volume 100 μL). The tubes were vortexed again, mixed gently for 24 hours at 32 °C and the reactions stopped by freezing at -80 °C. Samples of this reaction were analyzed as described above by SDS-PAGE gels

5 (Figure 167), IEF gels (Figure 168) and MALDI-TOF MS.

MALDI was also performed on the PEGylated rFSH. During ionization, SA-PEG is eliminated from the N-glycan structure of the glycoprotein. Native FSH gave a peak at 13928; AS-rFSH (13282); resialylated r-FSH (13332); PEG1000-rFSH (13515; 14960 (1); 16455 (2); 17796 (3); 19321 (4)); and PEG 10000 (23560 (1); 34790 (2); 45670 (3); and 10 56760 (4)).

## 25. Pharmacokinetic Study of GlycoPEGylated FSH

This example sets forth the *in vivo* testing of the pharmacokinetic properties glycoPEGylated Follicle Stimulating Hormone (FSH) prepared according to the methods of 15 the invention as compared to non-PEGylated FSH.

FSH, FSH-SA-PEG (1 kDa) and FSH-SA-PEG (10 kDa) were radioiodinated using standard conditions (Amersham Biosciences, Arlington Heights, IL) and formulated in phosphate buffered saline containing 0.1% BSA. After dilution in phosphate buffer to the appropriate concentration, each of the test FSH proteins (0.4 μg, each) was injected 20 intravenously into female Sprague Dawley rats (250-300 g body weight) and blood drawn at time points from 0 to 80 hours. Radioactivity in blood samples was analyzed using a gamma counter and the pharmacokinetics analyzed using standard methods (Figure 169). FSH was cleared from the blood much more quickly than FSH-PEG(1 kDa), which in turn was clear somewhat more quickly than FSH-PEG(10 kDa).

25

## 26. Sertoli Cell Bioassay for *In Vitro* Activity of GlycoPEGylated FSH

This example sets forth a bioassay for follicle stimulating hormone (FSH) activity based on cultured Sertoli cells. This assay is useful to determine the bioactivity of FSH after glycan remodeling, including glycoconjugation.

30 This bioassay is based on the dose-response relationship that exists between the amount of estradiol produced when FSH, but not lutenizing hormone (LH), is added to

cultured Sertoli cells obtained from immature old rats. Exogenous testosterone is converted to 17 $\beta$ -estradiol in the presence of FSH.

Seven to 10 days old Sprague-Dawley rats were used to obtain Sertoli cells. After sacrifice, testes were decapsulated and tissue was dispersed by incubation in collagenase (1 mg/ml), trypsin (1mg/ml), hyaluronidase (1 mg/ml) and DNases (5  $\mu$ g/ml) for 5 to 10 min. The tubule fragments settled to the bottom of the flask and were washed in PBS (1x). The tubule fragments were reincubated for 20 min with a media containing the same enzymes: collagenase (1 mg/ml), trypsin (1mg/ml), hyaluronidase (1 mg/ml) and DNases (5  $\mu$ g/ml).

The tubule fragments were homogenized and plated into a 24 well plate in a serum free media. 5 x 10<sup>5</sup> cells were dispersed per well. After 48h incubation at 37° C and 5% CO<sub>2</sub>, fresh media was added to the cells. Composition of the serum free media: DMEM (1 vol), Ham's F10 nutrient mixture (1 vol), insulin 1  $\mu$ g/ml, Transferrin 5  $\mu$ g/ml, EGF 10 ng/ml, T4 20 pg/ml, Hydrocortisone 10<sup>-8</sup> M, Retinoic acid 10<sup>-6</sup> M.

The stimulation experiment consists of a 24 hour incubation with standard FSH or samples at 37°C and 5% CO<sub>2</sub>. The mean intra-assay coefficient of variation is 9% and the mean inter-assay coefficient of variation is 11%.

The 17B-estradiol Elisa Kit DE2000 (R&D Systems, Minneapolis, MN) was used to quantify the level of estradiol after incubation with FSH, FSH-SA-PEG (1 kDa) and FSH-SA-PEG (10 kDa).

The procedure was as follows: 100  $\mu$ l of Estradiol Standard (provided with kit and prepared as per instructions with kit) or sample was pipetted into wells of 17B-estradiol Elisa plate(s); 50  $\mu$ l of 17B-estradiol Conjugate (provided with kit, prepared as per instructions with kit) was added to each well; 50  $\mu$ l of 17B-estradiol antibody solution (provided with kit and prepared as per instructions with kit) was added to each well; plates were incubated for 2 hour at room temperature at 200 rpm; the liquid was aspirated from each well; the wells were washed 4 times using the washing solution; all the liquid was removed from the wells; 200  $\mu$ l of pNPP Substrate (provided with kit and prepared as per instructions with kit) was added to all wells and incubated for 45 min; 50  $\mu$ l of Stop solution (provided with kit and prepared as per instructions with kit) was added and the plates were read it at 405 nm (Figure 170).

While FSH-PEG(10 kDa) exhibited a modest stimulation of Sertoli cells, at 1  $\mu$ g/ml, FSH-PEG(1 kDa) stimulated Sertoli cells up to 50% more than unPEGylated FSH.



## 27. Steelman-Pohley Bioassay of *In Vivo* Activity of GlycoPEGylated FSH

In this example, the Steelman-Pohley bioassay (Steelman and Pohley, 1953, Endocrinology 53:604-615) was used to determine the *in vivo* activity of glycoPEGylated FSH. The Steelman-Pohley assay uses the change in ovary weight of a rat to measure the *in vivo* activity of FSH that is coinjected with human chorionic gonadotropin.

The Steelman-Pohley bioassay was performed according to the protocol described in Christin-Maitre et al. (2000, Methods 21:51-57). Seventy female Sprague-Dawley Rats (Charles River Laboratories, Wilmington, MA), aged 21 to 22 days, were housed in the testing facility for at least 5 days before the beginning of the assay procedure. Throughout the procedure, the animal room was climate controlled at 18 to 26°C, 30 to 70% relative humidity, and 12 hr. artificial light/12 hr. dark. All animals were fed Certified Rodent Chow (Harlan Teklad, Madison WI) or the equivalent, and water, both *ad libitum*. Animal procedures were performed at Calvert Preclinical Services, Inc. (Olyphant, PA).

Recombinant FSH was expressed in CHO cells, purified by standard techniques and glycoPEGylated with PEG (1 kDa). The rats were divided into seven test groups, with ten animals per group. On days -1 and 0, animals of all groups were subcutaneously injected with 20 I.U. of human chorionic gonadotropin (HCG) in 0.5 ml of 0.9 % NaCl. On days 1, 2 and 3, the control animals were subcutaneously injected with a dose of 0.5 ml containing 20 I.U. HCG in 0.9% NaCl, while in the other groups, the HCG dose was augmented with either rFSH or rFSH-SA-PEG (1 kDa) at either 0.14 µg, 0.4 µg or 1.2 µg per dose. On day 4, the animals were euthanized by CO<sub>2</sub> inhalation. The ovaries were removed, trimmed and weighed. The average ovary weight was determined for each group.

Figure 171 presents the average ovary weight of the test groups on day 4. The groups receiving HCG alone (control) or the low dose (0.14 µg) of either rFSH or rFSH-SA-PEG (1 kDa) had ovary weights that were roughly equivalent. The groups receiving the medium (0.4 µg) or high (1.2 µg) doses of rFSH or rFSH-SA-PEG (1 kDa) had ovary weights roughly twice that of the control group. At the medium dose (0.4 µg), the glycoPEGylated rFSH had roughly the same *in vivo* activity (as determined by ovary weight) as the unPEGylated rFSH.

At the high dose (1.2 µg), the glycoPEGylated rFSH had somewhat higher *in vivo* activity than the unPEGylated rFSH.

### G-CSF

#### 28. GlycoPEGylation of G-CSF produced in CHO cells

**Preparation of Asialo-Granulocyte-Colony Stimulation Factor (G-CSF).** G-CSF produced in CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl<sub>2</sub> and concentrated to 500 µL in a Centricon Plus 20 centrifugal filter. The solution is incubated with 300 mU/mL Neuraminidase II (*Vibrio cholerae*) for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel performed. The reaction mixture is then added to prewashed N-(*p*-aminophenyl)oxamic acid-agarose conjugate (800 µL/mL reaction volume) and the washed beads gently rotated for 24 hours at 4 °C. The mixture is centrifuged at 10,000 rpm and the supernatant was collected. The beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer and once with 0.2 mL of the Tris-EDTA buffer and all supernatants are pooled. The supernatant is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub> and then twice more against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub>. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter and stored at -20 °C. The conditions for the IEF gel were run according to the procedures and reagents provided by Invitrogen. Samples of native and desialylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

**Preparation of G-CSF-(alpha2,3)-Sialyl-PEG.** Desialylated G-CSF was dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST3Gal1 at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis

according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

**Preparation of G-CSF-(alpha2,8)-Sialyl-PEG.** G-CSF produced in CHO cells, which contains an alpha2,3-sialylated O-linked glycan, is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of CST-II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction has CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

**Preparation of G-CSF-(alpha2,6)-Sialyl-PEG.** G-CSF, containing only O-linked GalNAc, is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST6GalNAcI or II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction has CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

G-CSF produced in CHO cells was treated with Arthrobacter sialidase and was then purified by size exclusion on Superdex75 and was treated with ST3GalI or ST3Gal2 and then with CMP-SA-PEG 20Kda. The resulting molecule was purified by ion exchange and

gel filtration and analysis by SDS PAGE demonstrated that the PEGylation was complete. This is the first demonstration of glycoPEGylation of an O-linked glycan.

### Glucocerebrosidase

#### 5        29. Glucocerebrosidase-mannose-6-phosphate produced in CHO cells

This example sets forth the procedure to glycoconjugate mannose-6-phosphate to a peptide produced in CHO cells such as glucocerebrosidase.

- Preparation of asialo-glucoceramide.** Glucocerebrosidase produced in CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, and is
- 10 incubated with 300 mU/mL sialidase-agarose conjugate for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel and SDS-PAGE performed according to Invitrogen procedures. The mixture is centrifuged at 10,000 rpm and the supernatant is collected. The beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer, and once with 0.2 mL of the Tris-EDTA buffer.
- 15 All supernatants are pooled. The supernatant is dialyzed at 4 °C against 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub> and then twice more against 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub>. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed
- 20 against water and analyzed by MALDI-TOF MS.

#### **Preparation of Glucocerebrosidase-SA-linker-Mannose-6-phosphate (procedure**

- 1).** Asialo-glucocerebrosidase from above is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-linker-Man-6-phosphate and 0.1 U/mL of ST3Gal3 at 32°C for 2 days. To monitor the
- 25 incorporation of sialic acid-linker-Man-6-phosphate, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas TSK-Gel-3000 analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. When the reaction is complete, the reaction mixture is purified
- 30 using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using

SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

**Preparation of Glucocerebrosidase-SA-linker-Mannose-6-phosphate (procedure**

- 2). Glucocerebrosidase, produced in CHO but incompletely sialylated, is dissolved at 2.5  
5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated  
with 1 mM CMP-sialic acid-linker-Man-6-phosphate and 0.1 U/mL of ST3Gal3 at 32°C for 2  
days. To monitor the incorporation of sialic acid-linker-Man-6-phosphate, a small aliquot of  
the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the  
peptide is separated from the free label by gel filtration on a Toso Haas TSK-Gel-3000  
10 analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the  
peptide is quantitated using an in-line fluorescent detector. When the reaction is complete,  
the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using  
PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the  
reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and  
15 reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by  
MALDI-TOF MS.

30. Glucocerebrosidase-transferrin

- This example sets forth the procedures for the glycoconjugation of proteins, and in  
20 particular, transferrin is glycoconjugated to glucocerebrosidase. The GlcNAc-ASN structures  
are created on glucoceraminidase, and Transferrin-SA-Linker-Gal-UDP is conjugated to  
GNDF GlcNAc-ASN structures using galactosyltransferase.

- Preparation of GlcNAc-glucocerebrosidase (Cerezyme™).** Cerezyme™  
(glucocerebrosidase) produced in CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM  
25 Tris-HCl pH 7.4, 0.15 M NaCl, and is incubated with 300 mU/mL Endo-H-agarose conjugate  
for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with  
the appropriate buffer and a IEF gel and SDS-PAGE performed according to Invitrogen  
procedures. The mixture is centrifuged at 10,000 rpm and the supernatant is collected. The  
beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer and  
30 once with 0.2 mL of the Tris-EDTA buffer and all supernatants are pooled. The supernatant  
is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub> and then twice

more against 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub>. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF

5 MS.

**Preparation of Transferrin-SA-Linker-Gal-glucocerebrosidase.** Transferrin-SA-Linker-Gal-UDP from above is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl<sub>2</sub>, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 2.5 mg/mL GlcNAc-glucocerebrosidase and 0.1 U/mL of galactosyltransferase at 32°C for 2 days. To monitor the incorporation of glucocerebrosidase, the peptide is separated by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1) and the product detected by UV absorption. The reaction mixture is then purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

## GM-CSF

### 31. Generation and PEGylation of GlcNAc-ASN Structures: GM-CSF produced in *Saccharomyces*

This example sets forth the preparation of Tissue-type Activator with PEGylated GlcNAc-Asn structures.

Recombinant GM-CSF expressed in yeast is expected to contain 2 N-linked and 2 O-linked glycans. The N-linked glycans should be of the branched mannan type. This recombinant glycoprotein is treated with an endoglycosidase from the group consisting of endoglycosidase H, endoglycosidase-F1, endoglycosidase-F2, endoglycosidase-F3, endoglycosidase-M either alone or in combination with mannosidases I, II and III to generate GlcNAc nubs on the asparagine (Asn) residues on the peptide/protein backbone.

The GlcNAc-Asn structures on the peptide/protein backbone is then be modified with galactose or galactose-PEG using UDP-galactose or UDP-galactose-6-PEG, respectively, and a galactosyltransferase such as GalT1. In one case the galactose-PEG is the terminal residue.

In the second case the galactose is further modified with SA-PEG using a CMP-SA-PEG donor and a sialyltransferase such as ST3GalIII. In another embodiment the GlcNAc-Asn structures on the peptide/protein backbone can be galactosylated and sialylated as described above, and then further sialylated using CMP-SA-PEG and an  $\alpha$ 2,8-sialyltransferase such as the enzyme encoded by the *Campylobacter jejuni* cst-II gene.

### Herceptin™

#### 32. Glycoconjugation of mithramycin to Herceptin™

This example sets forth the procedures to glycoconjugate a small molecule, such as mithramycin to Fc region glycans of an antibody molecule produced in mammalian cells. Here, the antibody Herceptin™ is used, but one of skill in the art will appreciate that the method can be used with many other antibodies.

**Preparation of Herceptin™-Gal-linker-mithramycin.** Herceptin™ is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl<sub>2</sub>, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 1 mM UDP-galactose-linker-mithramycin and 0.1 U/mL of galactosyltransferase at 32°C for 2 days to introduce the mithramycin in the Fc region glycans. To monitor the incorporation of galactose, a small aliquot of the reaction has <sup>14</sup>C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

When the reaction is complete, the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The fractions containing product are combined, concentrated, buffer exchanged and then freeze-dried. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Interferon  $\alpha$  and Interferon  $\beta$ 33. GlycoPEGylation of Proteins expressed in Mammalian or Insect Systems:  
EPO, Interferon  $\alpha$  and Interferon  $\beta$ 

5 This example sets forth the preparation of PEGylated peptides that are expressed in mammalian and insect systems.

**Preparation of acceptor from mammalian expression systems.** The peptides to be glycoPEGylated using CMP-sialic acid PEG need to have glycans terminating in galactose. Most peptides from mammalian expression systems will have terminal sialic acid that first  
10 needs to be removed.

**Sialidase digestion.** The peptide is desialylated using a sialidase. A typical procedure involves incubating a 1 mg/mL solution of the peptide in Tris-buffered saline, pH 7.2, with 5 mM  $\text{CaCl}_2$  added, with 0.2 U/mL immobilized sialidase from *Vibrio cholera* (Calbiochem) at 32°C for 24 hours. Microbial growth can be halted either by sterile filtration  
15 or the inclusion of 0.02% sodium azide. The resin is then removed by centrifugation or filtration, and then washed to recover entrapped peptide. At this point, EDTA may be added to the solution to inhibit any sialidase that has leached from the resin.

**Preparation from insect expression systems.** EPO, interferon-alpha, and interferon-beta may also be expressed in non-mammalian systems such as yeast, plants, or  
20 insect cells. The peptides to be glycoPEGylated using CMP-sialic acid PEG need to have glycans terminating in galactose. The majority of the N-glycans on peptides expressed in insect cells, for example, are the trimannosyl core. These glycans are first built out to glycans terminating in galactose before they are acceptors for sialyltransferase.

**Building acceptor glycans from trimannosyl core.** Peptide (1 mg/mL) in Tris-buffered saline, pH 7.2, containing 5 mM  $\text{MnCl}_2$ , 5 mM UDP-glcNAc, 0.05 U/mL GLCNACT I, 0.05 U/mL GLCNACT II, is incubated at 32°C for 24 hours or until the reaction is substantially complete. Microbial growth can be halted either by sterile filtration  
25 or the inclusion of 0.02% sodium azide. After buffer exchange to remove UDP and other small molecules, UDP-galactose and  $\text{MnCl}_2$  are each added to 5 mM, galactosyltransferase is added to 0.05 U/mL, and is incubated at 32°C for 24H or until the reaction is substantially  
30



complete. Microbial growth can be halted either by sterile filtration or the inclusion of 0.02% sodium azide. The peptides are then ready for glycoPEGylation.

**Building O-linked glycans.** A similar strategy may be employed for interferon alpha to produce enzymatically the desired O-glycan Gal-GalNAc. If necessary, GalNAc linked to serine or threonine can be added to the peptide using appropriate peptide GalNAc transferases (e.g. GalNAc T1, GalNAc T2, T3, T4, etc. ) and UDP-GalNAc. Also, if needed, galactose can be added using galactosyltransferase and UDP-galactose.

**GlycoPEGylation using sialyltransferase.** The glycopeptides (1 mg/mL) bearing terminal galactose in Tris buffered saline + 0.02% sodium azide are incubated with CMP-SA-PEG (0.75 mM) and 0.4 U/mL sialyltransferase (ST3Gal3 or ST3Gal4 for N-glycans on EPO and interferon beta; ST3Gal4, or ST3Gal1 for O-glycans on interferon alpha) at 32°C for 24 hours. Other transferases that may work include the 2,6 sialyltransferase from *Photobacterium damsella*. The acceptor peptide concentration is most preferably in the range of 0.1 mg/mL up to the solubility limit of the peptide. The concentration of CMP-SA-PEG should be sufficient for there to be excess over the available sites, but not so high as to cause peptide solubility problems due to the PEG, and may range from 50 µM up to 5 mM, and the temperature may range from 2°C up to 40°C. The time required for complete reaction will depend on the temperature, the relative amounts of enzyme to acceptor substrate, the donor substrate concentration, and the pH.

#### 34. GlycoPEGylation of Interferon $\alpha$ produced in CHO cells

**Preparation of Asialo-Interferon  $\alpha$ .** Interferon alpha produced from CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl<sub>2</sub> and concentrated to 500 µL in a Centricon Plus 20 centrifugal filter. The solution is incubated with 300 mU/mL Neuraminidase II (*Vibrio cholerae*) for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel performed. The reaction mixture is then added to prewashed N-(*p*-aminophenyl)oxamic acid-agarose conjugate (800 µL/mL reaction volume) and the washed beads gently rotated for 24 hours at 4 °C. The mixture is centrifuged at 10,000 rpm and the supernatant was collected. The beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer and once with 0.2 mL of the Tris-EDTA buffer and all

supernatants were pooled. The supernatant is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub> and then twice more against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub>. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter and stored at -20 °C. The conditions for the IEF gel are run according to the procedures and reagents provided by Invitrogen. Samples of native and desialylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

**Preparation of Interferon-alpha-(alpha2,3)-Sialyl-PEG.** Desialylated interferon-alpha is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST3Gal1 at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and desialylated Interferon-alpha are dialyzed against water and analyzed by MALDI-TOF MS.

**Preparation of Interferon-alpha-(alpha2,8)-Sialyl-PEG.** Interferon-alpha produced in CHO, which contains an alpha2,3-sialylated O-linked glycan, is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of CST-II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction has CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis

according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated interferon-alpha are dialyzed against water and analyzed by MALDI-TOF MS.

**Preparation of Interferon-alpha-(alpha2,6)-Sialyl-PEG.** Interferon-alpha, containing only O-linked GalNAc, was dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST6GalNAcI or II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated interferon-alpha are dialyzed against water and analyzed by MALDI-TOF MS.

### 35. GlycoPEGylation of Interferon-β-1a with PEG (10 kDa) and PEG (20 kDa)

This example illustrates a procedure PEGylate Interferon-β with either PEG (10 kDa) or PEG (20 kDa).

Briefly, Interferon-β-1a (INF-β) was obtained from Biogen (Avonex™). The INF-β was first purified by Superdex-75 chromatography. The INF-β was then desialylated with *Vibrio cholerae* sialidase. The INF-β was then PEGylated with SA-PEG (10 kDa) or SA-PEG (20 kDa) and purified with Superdex-200 chromatography.

**Superdex-75 chromatography purification.** INF-β (150 μg) was applied to a Superdex-75 column (Amersham Biosciences, Arlington Heights, IL) and eluted with PBS with 0.5 M NaCl, 0.02 Tween-20, 20 mM histidine and 10% glycerol. The eluant was monitored for absorbance at 280 nm (Figure 172A and 172B) and fractions were collected. Peaks 4 and 5 were pooled, concentrated in an Amicon Ultra 15 spin filter (Millipore, Billerica, MA), and the buffer was exchanged to TBS with 5 mM CaCl<sub>2</sub>, 0.02% Tween-20, 20 mM histidine and 10% glycerol.

**Sialidase Reaction.** The INF- $\beta$  was then desialylated with *Vibrio cholera* sialidase (70 mU/ml, CALBIOCHEM®, EMD Biosciences, Inc., San Diego, CA) on agarose in TBS with 5 mM CaCl<sub>2</sub>, 0.02% Tween-20, 20 mM histidine and 10% glycerol. The reaction was carried out at 32°C for 18 hours. The INF- $\beta$  was removed from the agarose with a 0.22  $\mu$ m Spin-X™ filter (Corning Technology, Inc., Norcross, GA). Figure 173A depicts the MALDI analysis of glycans released from native INF- $\beta$ . The native INF- $\beta$  has many glycoforms containing terminal sialic acid moieties. Figure 173B depicts the MALDI analysis of glycans released from desialylated INF- $\beta$ . The desialylated INF- $\beta$  has primarily one glycoform which is bi-antennary with terminal galactose moieties.

**Lectin Dot-Blot Analysis of Sialylation.** Samples of the INF- $\beta$  from the desialidase reaction were dot-blotted onto nitrocellulose and then blocked with Tris buffered saline (TBS: 0.05M Tris, 0.15M NaCl, pH 7.5) and DIG kit (glycan differentiation kit available from Roche #1 210 238) blocking buffer. Some of the blots were incubated with *Maackia amurensis* agglutinin (MAA) labeled with digoxigenin (DIG) (Roche Applied Science, Indianapolis, IL) to detect  $\alpha$ 2,3-sialylation of INF- $\beta$ . These blots were washed with TBS then incubated with anti-digoxigenin antibody labeled with alkaline phosphatase, then washed again with TBS and developed with NBT/X-phosphate solution, wherein NBT is 4-nitro blue tetrazolium chloride and X-phosphate is 5-bromo-4-chloro-3-indoyl phosphate. The left side of Figure 174 depicts the results of the MAA blot of INF- $\beta$  after the desialylation reaction. The INF- $\beta$  is partially desialylated, as indicated by the decrease in dot development as compared to native INF- $\beta$  in the desialylated samples.

Other blots were incubated with *Erthrina cristagalli* lectin (ECL) labeled with biotin (Vector Laboratories, Burlingame, CA) to detect exposed galactose residues on INF- $\beta$ . After incubation with 2.5  $\mu$ g/ml ECL, the blots were washed in TBS and incubated with streptavidin labeled with alkaline phosphatase. The blots were then washed again and developed. The right side of Figure 174 depicts the ECL blot after development. The increased intensity of the dot of desialylated INF- $\beta$  as compared to the native INF- $\beta$  indicate more exposed galactose moieties and therefore extensive desialylation.

**PEGylation of Desialylated INF- $\beta$  with SA-PEG (10 kDa).** Desialylated INF- $\beta$  (0.05 mg/ml) was PEGylated with ST3Gal3 (50 mU/ml) and CMP-SA-PEG (10 kDa) (250

μM) in an appropriate buffer of TBS + 5 mM CaCl<sub>2</sub>, 0.02% Tween 20, 20 mM histidine, 10% glycerol for 50 hours at 32°C. Figure 175 depicts the SDS-PAGE analysis of the reaction products showing PEGylated INF-β at approximately 98 kDa.

**PEGylation of Desialylated INF-β with SA-PEG (20 kDa).** Desialylated INF-β (0.5 mg/ml) was PEGylated with ST3Gal3 (170 mU/ml) and CMP-SA-PEG (20 kDa) in an appropriate buffer of TBS + 5 mM CaCl<sub>2</sub>, 0.02% Tween 20, 20 mM histidine, 10% glycerol for 50 hours at 32°C. Figure 176 depicts the SDS-PAGE analysis the products of the PEGylation reaction. The PEGylated INF-β has many higher molecular weight bands not found in the unmodified INF-β indicating extensive PEGylation.

**Superdex-200 Purification of INF-β PEGylated with PEG (10 kDa).** The products of the PEGylation reaction were separated on a Superdex-200 column (Amersham Biosciences, Arlington Heights, IL) in PBS with 0.5 NaCl, 0.02 Tween-20, 20 mM histidine and 10% glycerol at 1ml/min and 30 cm/hr flow. The eluant was monitored for absorbance at 280 nm (Figure 177) and fractions were collected. Peaks 3 and 4 were pooled and concentrated in an Amicon Ultra 15 spin filter.

**Bioassay of INF-β PEGylated with PEG (10 kDa).**

The test is inhibition of the proliferation of the lung carcinoma cell line, A549. The A549 cell line are lung carcinoma adherent cells growing in RPMI + 10% FBS at 37°C 5% CO<sub>2</sub>. They can be obtained from ATCC # CCL-185. Wash the cells with 10 ml of PBS and remove the PBS. Add 5 ml of trypsin, incubate for 5 minutes at room temperature or 2 minutes at 37°C. When the cells are detached resuspend into 25 ml of media and count the cells. Dilute the cells at a concentration of 10000 cells/ml and add 200 ul / well (96 wells plate). Incubate for 4 hours at 37°C 5% CO<sub>2</sub>. Prepare 1 ml of IFN B at a concentration of 0.1 ug/ml. Filter it under the hood with a 0.2 um filter. Add 100 ul per well (8 replicates = 1 lane). Incubate for 3 days (do not let the cells go to confluence). Remove 200 ul of media (only 100ul per well left). Add 25 μl of MTT (Sigma) (5 mg/ml filtered 0.22μm). Incubate for 4 hours at 37°C and 5% CO<sub>2</sub>. Aspirate the media gently and add 100 μl of a mixture of isopropanol (100 ml and 6N HCl. Aspirate up and down to homogenize the crystal violet. Read OD 570nm (remove the background at 630 or 690 nm).

Figure 178 depicts the results of the bioassay of the peaks containing INF- $\beta$  PEGylated with PEG (10 kDa) as eluted from the Superdex-200 column.

**Superdex-200 Purification of INF- $\beta$  PEGylated with PEG (20 kDa).** The products of the PEG (20 kDa) PEGylation reaction were separated on a Superdex-200 column (Amersham Biosciences, Arlington Heights, IL) in PBS with 0.5 NaCl, 0.02 Tween-20, 20 mM histidine and 10% glycerol at 1 ml/min flow. The eluant was monitored for absorbance at 280 nm (Figure 179) and fractions were collected. Peak 3 contained most of the INF- $\beta$  PEGylated with PEG (20 kDa).

**Endotoxin test of INF- $\beta$  PEGylated with PEG (20 kDa).**

Limulus Lysate Test was performed, BioWhittaker # 50-647U

**Table 24.** Results of the endotoxin test of INF- $\beta$  PEGylated with PEG (20 kDa).

	Concentration		
INF- $\beta$ with PEG (20 kDa)	10 EU/ml	0.06 mg/ml	0.16 EU/ $\mu$ g
INF- $\beta$ with PEG (20 kDa)	1 EU/ml	0.07 mg/ml	0.014 EU/ $\mu$ g
Native INF- $\beta$	40 EU/ml	0.1 mg/ml	0.4 EU/ $\mu$ g

**Remicade™**

**36. GlycoPEGylation of Remicade™ antibody**

This example sets forth the procedure to glycoPEGylate a recombinant antibody molecule by introducing PEG molecules to the Fc region glycans. Here Remicade™, a TNF-R: IgG Fc region fusion protein, is the exemplary peptide.

**Preparation of Remicade™-Gal-PEG (10 kDa).** Remicade™ is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl<sub>2</sub>, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 1 mM UDP-galactose-PEG (10 kDa) and 0.1 U/mL of galactosyltransferase at 32°C for 2 days to introduce the PEG in the Fc region glycans. To monitor the incorporation of galactose, a small aliquot of the reaction has <sup>14</sup>C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

When the reaction is complete, the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The fractions containing product are combined, concentrated, buffer exchanged and then freeze-dried. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

### Rituxan™

#### 37. Glycoconjugation of geldanamycin to Rituxan™

This example sets forth the glycoconjugation of a small molecule, such as geldanamycin, to the Fc region glycans of an antibody produced in CHO cells, such as Rituxan™. Here, the antibody Rituxan™ is used, but one of skill in the art will appreciate that the method can be used with many other antibodies.

**Preparation of Rituxan™-Gal-linker-geldanamycin.** Rituxan™ is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl<sub>2</sub>, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 1 mM UDP-galactose-linker-geldanamycin and 0.1 U/mL of galactosyltransferase at 32°C for 2 days to introduce the geldanamycin in the Fc region glycans. To monitor the incorporation of galactose, a small aliquot of the reaction has <sup>14</sup>C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

When the reaction is complete, the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The fractions containing product are combined, concentrated, buffer exchanged and then freeze-dried. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Rnase38. Remodeling high mannose N-glycans to hybrid and complex N-glycans:  
Bovine pancreatic RNase

This example sets forth the preparation of bovine pancreas RNase with hybrid or  
5 complex N-glycans. The high mannose N-linked glycans of the RNase are enzymatically  
digested and elaborated to create hybrid N-linked glycans. Additionally, the high mannose  
N-linked glycans of the RNase are enzymatically digested and elaborated to create complex  
N-linked glycans.

High mannose structures of *N*-linked oligosaccharides in glycopeptides can be  
10 modified to hybrid or complex forms using the combination of  $\alpha$ -mannosidases and  
glycosyltransferases. This example summarizes the results in such efforts using a simple *N*-  
Glycan as a model substrate.

Ribonuclease B (RNaseB) purified from bovine pancreas (Sigma) is a glycopeptide  
consisting of 124 amino acid residues. It has a single potential *N*-glycosylation site modified  
15 with high mannose structures. Due to its simplicity and low molecular weight (13.7 kDa to  
15.5 kDa), ribonuclease B is a good candidate to demonstrate the feasibility of the *N*-Glycan  
remodeling from high mannose structures to hybrid or complex *N*-linked oligosaccharides.  
The MALDI-TOF spectrum of RNaseB (Figure 180A) and HPLC profile for the  
oligosaccharides cleaved from RNaseB by N-Glycanase (Figure 180B) indicated that, other  
20 than a small portion of the non-modified peptide, the majority of *N*-glycosylation sites of the  
peptide are modified with high mannose oligosaccharides consisting of 5 to 9 mannose  
residues.

**Conversion of high mannose N-Glycans to hybrid N-Glycans.** High mannose *N*-  
Glycans were converted to hybrid *N*-Glycans using the combination of  $\alpha$ 1,2-mannosidase,  
25 GlcNAcT-I ( $\beta$ -1,2-*N*-acetyl glucosaminyl transferase), GalT-I ( $\beta$ 1,4-galactosyltransferase) and  
 $\alpha$ 2,3-sialyltransferase /or  $\alpha$ 2,6-sialyltransferase as shown in Figure 181.

As an example, high mannose structures in RNaseB were successfully converted to  
hybrid structures.

Man<sub>5</sub>GlcNAc<sub>2</sub>-R was obtained from Man<sub>5-9</sub>GlcNAc<sub>2</sub>-R catalyzed by a single  $\alpha$ 1,2-  
30 mannosidase cloned from *Trichoderma reesei* (Figure 182). RNase B (1 g, about 67  $\mu$ mol)  
was incubated at 30°C for 45 hr with 15 mU of the recombinant *T. reesei*  $\alpha$ 1,2-mannosidase



in MES buffer (50 mM, pH 6.5) in a total volume of 10 mL.  $\text{Man}_6\text{GlcNAc}_2$ -protein structures have been successfully converted to  $\text{Man}_5\text{GlcNAc}_2$ -protein with high efficiency by the recombinant mannosidase.

Alternately,  $\text{Man}_5\text{GlcNAc}_2$ -R was obtained from  $\text{Man}_5\text{GlcNAc}_2$ -R catalyzed by a single  $\alpha$ 1,2-mannosidase purified from *Aspergillus saitoi* (Figure 183). RNase B (40  $\mu$ g, about 2.7 nmol) was incubated at 37°C for 42.5 hr with 25  $\mu$ U of the commercial *A. saitoi*  $\alpha$ 1,2-mannosidase (Glyko or CalBioChem) in NaOAc buffer (100 mM, pH 5.0) in a total volume of 20  $\mu$ L.  $\text{Man}_6\text{GlcNAc}_2$ -protein structures were successfully converted to  $\text{Man}_5\text{GlcNAc}_2$ -protein by the commercially available mannosidase. However, a new peak corresponding to the GlcNAc-protein appears in the spectrum, indicating the possible contamination of endoglycosidase H in the preparation. Although several mammalian alpha-mannosidases were required to achieve this step, the fungal  $\alpha$ 1,2-mannosidase was very efficient to remove all  $\alpha$ 1,2-linked mannose residues.

GlcNAcT-I then added a GlcNAc residue to the  $\text{Man}_5\text{GlcNAc}_2$ -R (Figure 184). The reaction mixture after the *T. reesei*  $\alpha$ 1,2-mannosidase reaction containing RNase B (600  $\mu$ g, about 40 nmol) was incubated with non-purified recombinant GlcNAcT-I (34 mU) in MES buffer (50 mM, pH 6.5) containing  $\text{MnCl}_2$  (20 mM) and UDP-GlcNAc (5 mM) in a total volume of 400  $\mu$ L at 37°C for 42 hr. A GlcNAc residue was quantitatively added to  $\text{Man}_5\text{GlcNAc}_2$ -protein by the recombinant GlcNAcT-I.

A Gal residue was then added using GalT 1 (Figure 185). The reaction mixture after the GnT-I reaction containing RNase B (120  $\mu$ g, about 8 nmol) was incubated at 37°C for 20 hr with 3.3 mU of the recombinant GalT-1 in Tris-HCl buffer (100 mM, pH 7.3) containing UDP-Gal (7.5 mM) and  $\text{MnCl}_2$  (20 mM) in a total volume of 100  $\mu$ L. A Gal residue was added to about 98% of the GlcNAc- $\text{Man}_5\text{GlcNAc}_2$ -protein by the recombinant GalT 1.

The next step was the addition of a sialic acid using an  $\alpha$ 2,3-sialyltransferase or an  $\alpha$ 2,6-sialyltransferase (Figure 186). As an example, ST3Gal III, an  $\alpha$ 2,3-sialyltransferase was used. The reaction mixture after the GalT-1 reaction containing RNase B (13  $\mu$ g, about 0.87 nmol) was incubated at 37°C for 16 hr with 8.9 mU of recombinant ST3Gal III in Tris-HCl buffer (100 mM, pH 7.3) containing CMP-Sialic acid (5 mM) and  $\text{MnCl}_2$  (20 mM) in a total volume of 20  $\mu$ L. A sialic acid residue was added to about 90% of the Gal-GlcNAc-

Man<sub>5</sub>GlcNAc<sub>2</sub>-protein by recombinant ST3Gal III using CMP-SA as the donor. The yield can be further improved by adjusting the reaction conditions.

For convenience, no purification or dialysis step was required after each reaction described above. More interesting, GalT I and ST3Gal III can be combined in a one-pot reaction. Similar yields were obtained as compared with the separate reactions. The reaction mixture after the GlcNAcT-I reaction containing RNase B (60 µg, about 4 nmol) was incubated at 37°C for 20 hr with 1.7 mU of recombinant GalT I, 9.8 mU of recombinant ST3Gal III in Tris-HCl buffer (100 mM, pH 7.3) containing UDP-Gal (7.5 mM), CMP-sialic acid (5 mM) and MnCl<sub>2</sub> (20 mM) in a total volume of 60 µl.

As shown in Figure 187, SA-PEG (10 kDa) was successfully added to the RNaseB. The reaction mixture after the GalT-I reaction containing RNase B (6.7 µg, about 0.45 nmol) was dialyzed against H<sub>2</sub>O for 1 hour at room temperature and incubated at 37°C for 15.5 hours with 55 mU of the recombinant ST3Gal III in Tris-HCl buffer (50 mM, pH 7.3) containing CMP-SA-PEG (10 kDa) (0.25 mM) and MnCl<sub>2</sub> (20 mM) in a total volume of 20 µl. PEG-modified sialic acid residues were successfully added to the Gal-GlcNAc-Man<sub>5</sub>GlcNAc<sub>2</sub>-peptide by the recombinant ST3Gal III. The yield can be further improved by adjusting the reaction conditions.

**Conversion of high mannose N-Glycans to complex N-Glycans.** To achieve this conversion, a GlcNAcβ1,2Man<sub>5</sub>GlcNAc<sub>2</sub>-peptide intermediate is obtained. As shown in Figure 188, there are at least four feasible routes to carry out the reaction from Man<sub>5</sub>GlcNAc<sub>2</sub>-peptide to this intermediate:

**Route I:** The Man<sub>5</sub>GlcNAc<sub>2</sub>-peptide produced by the fungal α1,2 mannosidase is a substrate of GlcNAc transferase I (GlcNAcT-I, enzyme 2) which adds one GlcNAc. The terminal α1,3- and α1,6-linked mannose residues of GlcNAcMan<sub>5</sub>GlcNAc<sub>2</sub>-peptide is removed by Golgi α-mannosidase II (ManII, enzyme 5). This route is a part of the natural pathway for the processing of N-linked oligosaccharides carried out in higher organisms.

**Route II:** Two mannose residues are first removed by an α-mannosidase (enzyme 6), then a GlcNAc is added by GlcNAcT-I (enzyme 2). Other than its natural acceptor Man<sub>5</sub>GlcNAc<sub>2</sub>-R, GlcNAcT-I can also recognize Man<sub>3</sub>GlcNAc<sub>2</sub>-R as its substrate and add one GlcNAc to the mannose core structure to form GlcNAcMan<sub>3</sub>GlcNAc<sub>2</sub>-peptide.

**Route III:** The  $\alpha$ 1,6-linked mannose is removed by an  $\alpha$ 1,6-mannosidase, followed by the addition of GlcNAc by GlcNAcT-I and removal of the terminal  $\alpha$ 1,3-linked mannose by an  $\alpha$ 1,3-mannosidase. From the experimental data obtained, GlcNAcT-I can recognize this Man<sub>4</sub>GlcNAc<sub>2</sub>-peptide as acceptor and add one GlcNAc residue to form

5 GlcNAcMan<sub>4</sub>GlcNAc<sub>2</sub>-peptide.

**Route IV:** Similar to Route III,  $\alpha$ 1,3-linked mannose is removed by an  $\alpha$ 1,3-mannosidase, followed by GlcNAcT-I reaction. Then the terminal  $\alpha$ 1,6-linked mannose can be removed by an  $\alpha$ 1,6-mannosidase.

After the function of GlcNAcT-I (responsible for the addition of the GlcNAc  $\beta$ 1,2-linked to the  $\alpha$ 1,3-mannose on the mannose core) and GlcNAcT-II (responsible for the addition of a second GlcNAc  $\beta$ 1,2-linked to the  $\alpha$ 1,6-mannose on the mannose core), the GlcNAc<sub>2</sub>Man<sub>3</sub>GlcNAc<sub>2</sub>-peptide can be processed by GalT 1 and sialyltransferase to form bi-antennary complex N- Glycans. Other GlcNAc transferases such as GlcNAcT-IV, GlcNAcT-V, and/or GlcNAcT-VI (Figure 188 and Figure 189) can also glycosylate the

10 linked to the  $\alpha$ 1,3-mannose on the mannose core) and GlcNAcT-II (responsible for the addition of a second GlcNAc  $\beta$ 1,2-linked to the  $\alpha$ 1,6-mannose on the mannose core), the GlcNAc<sub>2</sub>Man<sub>3</sub>GlcNAc<sub>2</sub>-peptide can be processed by GalT 1 and sialyltransferase to form bi-antennary complex N- Glycans. Other GlcNAc transferases such as GlcNAcT-IV, GlcNAcT-V, and/or GlcNAcT-VI (Figure 188 and Figure 189) can also glycosylate the

15 GlcNAc<sub>2</sub>Man<sub>3</sub>GlcNAc<sub>2</sub>-peptide. Additional glycosylation by the GalT 1 and sialyltransferases will form multi-antennary complex N-glycans. The enzyme GlcNAcT-III catalyzes the insertion of a bisecting GlcNAc, thus preventing the actions of ManII and subsequent action of transferases GlcNAcT-II, GlcNAcT-IV and GlcNAcT-V.

## 20 Tissue-Type Plasminogen Activator (TPA)

### 39. Fucosylation of TPA to create Sialyl Lewis X

This example sets forth the preparation of Tissue Tissue-type Plasminogen Activator (TPA) with N-linked sialyl Lewis X antigen.

**Sialylation.** TPA expressed in mammalian cells will often contain a majority of the glycans terminating in sialic acid, but to ensure complete sialylation, it would be beneficial to

25 first perform an *in vitro* sialylation. TPA in a suitable buffer (most preferably between pH 5.5 and 9, for example Tris buffered saline, pH 7.2) is incubated with CMP sialic acid and sialyltransferase for a time sufficient to convert any glycans lacking sialic acid to sialylated species. Typical conditions would be 1 mg/mL TPA, 3 mM CMP sialic acid, 0.02 U/mL

30 ST3Gal3, 32°C for 24 hours. Microbial growth can be halted either by sterile filtration or the

inclusion of 0.02% sodium azide. The TPA concentration is most preferably in the range 0.1 mg/mL up to the solubility limit of the peptide. The concentration of CMP-SA should be sufficient for there to be excess over the available sites, and might range from 50  $\mu$ M up to 50 mM, and the temperature from 2°C up to 40°C. The time required for complete reaction will depend on the temperature, the relative amounts of enzyme to acceptor substrate, the donor substrate concentration, and the pH. Other sialyltransferases that may be capable of adding sialic acid in 2,3 linkage include ST3Gal4; microbial transferases could also be used.

**Fucosylation.** Typical conditions for fucosylation would be 1 mg/mL TPA, 3 mM GDP-fucose, 0.02 U/mL FTVI, 5 mM MnCl<sub>2</sub>, 32°C for 24H in Tris buffered saline.

Microbial growth can be halted either by sterile filtration or the inclusion of 0.02% sodium azide. The TPA concentration is most preferably in the range 0.1 mg/mL up to the solubility limit of the peptide. The concentration of GDP-fucose should be sufficient for there to be excess over the available sites, and might range from 50  $\mu$ M up to 50 mM, and the temperature from 2°C up to 40°C. The time required for complete reaction will depend on the temperature, the relative amounts of enzyme to acceptor substrate, the donor substrate concentration, and the pH. Other fucosyltransferases that may be capable of making sialyl Lewis x include FTVII, FTV, FTIII, as well as microbial transferases could also be used.

#### 40. Trimming of high mannose to tri-mannose core structure: Tissue-type Plasminogen Activator produced in CHO

This example sets forth the preparation of Tissue-type Plasminogen Activator with a trimannose core by trimming back from a high mannose glycan.

Tissue-type plasminogen activator (TPA) is currently produced in Chinese Hamster Ovary (CHO) cells and contains a low amount of high mannose N-linked oligosaccharide.

The mannoses can be trimmed down using a variety of the specific mannosidases. The first step is to generate Man5GlcNAc2(Fuc0-1) from Man9GlcNAc2(Fuc0-1). This can be done using mannosidase I. Then either GlcNAcT1 (GlcNAc transferase I) is used to make GlcNAc1Man5GlcNAc2(Fuc0-1) or Mannosidase III is used to make Man3GlcNAc2(Fuc0-1). From Man3GlcNAc2(Fuc0-1), GlcNAc1Man3GlcNAc2(Fuc0-1) can be produced using GlcNAcT1 or from GlcNAc1Man5GlcNAc2(Fuc0-1), GlcNAc1Man3GlcNAc2(Fuc0-1) can be produced using Mannosidase II. GlcNAc1Man3GlcNAc2(Fuc0-1) is then converted into

GlcNAc2Man3GlcNAc2(Fuc0-1) using GlcNAcTransferase II (GlcNAcTII). The two terminal GlcNAc residues are then galactosylated using GalTI and then sialylated with SA-PEG using ST3GalIII.

- Conversely, TPA can be produce in yeast or fungal systems. Similar processing  
 5 would be required for fungal derived material.

#### 41. Generation and PEGylation of GlcNAc-ASN structures: TPA produced in Yeast

- This example sets forth the preparation of PEGylated GlcNAc-Asn structures on a  
 10 peptide such as TPA expressed in yeast.

Yeast expression is expected to result in a TPA which contains a single N-linked mannan-type structure. This recombinant glycoprotein is first treated with endoglycosidase H to generate GlcNAc structures on the asparagine (Asn) residues on the peptide.

- The GlcNAc-Asn structures on the peptide/protein backbone are then be modified  
 15 with galactose or galactose-PEG using UDP-galactose or UDP-galactose-6-PEG, respectively, and a galactosyltransferase such as GalT1. In one case, the galactose-PEG is the terminal residue. In the second case, the galactose is further modified with SA-PEG using a CMP-SA-PEG donor and a sialyltransferase such as ST3GalIII. In another embodiment, the GlcNAc-Asn structures on the peptide/protein backbone may be  
 20 galactosylated and sialylated as described above, and then further sialylated using CMP-SA-PEG and an  $\alpha$ 2,8-sialyltransferase such as the enzyme encoded by the *Campylobacter jejuni* cst-II gene.

#### Transferrin

- 25 42. GlycoPEGylation of Transferrin

This example sets forth the preparation of asialotransferrin and its sialylation with PEG-CMP-sialic acid.

- Preparation of Asialo-transferrin.** Human-derived holo-Transferrin, (10 mg) was dissolved in 500  $\mu$ L of 50 mM NaOAc, 5 mM CaCl<sub>2</sub>, pH 5.5. To this solution was added  
 30 500 mU Neuraminidase II (*Vibrio cholerae*) and the reaction mixture was shaken gently for 20.5 hours at 37 °C. The reaction mixture was added to the prewashed N-(p-

aminophenyl)oxamic acid-agarose conjugate (600  $\mu$ L) and the washed beads gently rotated for 24 hours at 4 °C. The mixture was centrifuged at 10,000 rpm and the supernatant was collected. The reaction mixture was adjusted to 5 mM EDTA by addition of 100  $\mu$ L of 30 mM EDTA to the washed beads, which were gently rotated for 20 hours at 4 °C. The suspension was centrifuged for 2 minutes at 10,000 rpm and the supernatant was collected. The beads were washed 5 times with 0.35 mL of 50 mM NaOAc, 5 mM CaCl<sub>2</sub>, 5 mM EDTA, pH 5.5 and all supernatants were pooled. The enzyme solution was dialyzed twice at 4 °C into 15 mM Tris-HCl, 1 M NaCl, pH 7.4. 0.3 mL of the transferrin solution (3.3 mL total) was removed and dialyzed twice against water. The remainder was dialyzed twice more at 4 °C against phosphate buffered saline. The dialyzed solution was stored at -20 °C. Protein samples were analyzed by IEF Electrophoresis. Samples (9  $\mu$ L, 25  $\mu$ g) were diluted with 16  $\mu$ L Tris buffer and mixed with 25  $\mu$ L of the sample loading buffer and applied to Isoelectric Focusing Gels (pH 3-7). Gels were run and fixed using standard procedures. Gels were stained with Colloidal Blue Stain.

**Sialyl-PEGylation of asialo-Transferrin.** Desialylated transferrin (250  $\mu$ g) and CMP-sialic acid or CMP-SA-PEG (1 kDa or 10 kDa)(0.05  $\mu$ mol) were dissolved in 69  $\mu$ L 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2 in 1.5 mL plastic tubes. The tubes were vortexed briefly and 100 mU ST3Gal3 (90  $\mu$ L) were added (total volume 250  $\mu$ L). The tubes were vortexed again and mixed gently for 24 hours at 32 °C. The reactions were stopped by freezing at -80 °C. Novex Tris-Glycine 8-16% 1 mm gels were used for SDS PAGE analysis (Figure 190). Samples (25  $\mu$ L, 25  $\mu$ g) were mixed with 25  $\mu$ L of sample loading buffer and 0.4  $\mu$ L of  $\beta$ -mercaptoethanol and heated for 6 minutes at 85 °C. Gels were run using standard conditions and stained with Colloidal Blue Stain. IEF gels were also performed as described above (Figure 191). Samples were also dialyzed against water and analyzed by MALDI-TOF.

**Results.** MALDI was also performed. Native transferrin (78729); asialotransferrin (78197); resialylated transferrin (79626/80703); with SA-PEG 1k (79037 (1); 80961 (2); 82535 (3); 84778 (4)); with SA-PEG 5k (90003 (2); 96117 (3); 96117 (4)); with SA-PEG 10k (100336 (2); 111421 (3); 122510 (4)).

#### 43. Transferrin-GDNF

This example sets forth the procedures for the glycoconjugation of proteins, and in particular, transferrin is glycoconjugated to GDNF. Transferrin-SA-Linker-Gal-UDP is prepared from transferrin. The galactose residue is removed from GDNF glycans, and Transferrin-SA-Linker-Gal-UDP is conjugated to GDNF glycans using a galactosyltransferase.

**Preparation of agalacto-GDNF.** GDNF produced in NSO cells (NSO murine myeloma cells) is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, and is incubated with 300 mU/mL beta-galactosidase-agarose conjugate for 16 hours at 32°C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel performed according to Invitrogen procedures. The mixture is centrifuged at 10,000 rpm and the supernatant is collected. The supernatant is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub> and then twice more against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN<sub>3</sub>. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter and stored at -20 °C. The conditions for the IEF gel are run according to the procedures and reagents provided by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

**Preparation of Transferrin-SA-Linker-Gal-UDP.** Asialo-transferrin is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with CMP-sialic acid-linker-Gal-UDP (molar amount to add 1 molar equivalent of nucleotide sugar to transferrin) and 0.1 U/mL of ST3Gal3 at 32°C for 2 days. To monitor the incorporation of sialic acid, a small aliquot of the reaction has <sup>14</sup>C-SA-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

The solution is incubated with 5 mM CMP-sialic acid and 0.1 U/mL of ST3Gal3 (to cap any unreacted transferrin glycans) at 32°C for 2 days. The incorporation into the peptide is quantitated using an in-line UV detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE

and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

- Preparation of Transferrin-SA-Linker-Gal-GDNF.** The transferrin-SA-Linker-Gal-UDP prepared as described above is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl<sub>2</sub>, 0.05% NaN<sub>3</sub>, pH 7.2. The solution is incubated with 2.5 mg/mL agalacto-GDNF and 0.1 U/mL of galactosyltransferase at 32°C for 2 days. To monitor the incorporation of galactose, a small aliquot of the reaction has <sup>14</sup>C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

- When the reaction is complete, the solution is incubated with 5 mM UDP-Gal and 0.1 U/mL of galactosyltransferase (to cap any unreacted transferrin glycans) at 32°C for 2 days followed by addition of 5 mM CMP-SA and 0.1 U/mL of ST3Gal3. After 2 additional days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

- The disclosures of each and every patent, patent application, and publication cited herein are hereby incorporated herein by reference in their entirety.

- While this invention has been disclosed with reference to specific embodiments, it is apparent that other embodiments and variations of this invention may be devised by others skilled in the art without departing from the true spirit and scope of the invention. The appended claims are intended to be construed to include all such embodiments and equivalent variations.



What is claimed:

1. An EPO peptide comprising one or more glycans, having a glycoconjugate molecule covalently attached to said peptide.

2. The EPO peptide of claim 1, wherein said one or more glycans is a  
5 monoantennary glycan.

3. The EPO peptide of claim 1, wherein said one or more glycans is a biantennary glycan.

4. The EPO peptide of claim 1, wherein said one or more glycans is a triantennary glycan.

10 5. The EPO peptide of claim 1, wherein said one or more glycans is at least a triantennary glycan.

6. The EPO peptide of claim 1, wherein said one or more glycans comprises at least two glycans comprising a mixture of mono or multiantennary glycans.

15 7. The EPO peptide of claim 1, wherein said one or more glycans is selected from an N-linked glycan and an O-linked glycan.

8. The EPO peptide of claim 1, wherein said one or more glycans is at least two glycans selected from an N-linked and an O-linked glycan.

9. The EPO peptide of claim 1, wherein said peptide is expressed in a cell selected from the group consisting of a prokaryotic cell and a eukaryotic cell.

20 10. The EPO peptide of claim 9, wherein said eukaryotic cell is selected from the group consisting of a mammalian cell, an insect cell and a fungal cell.

11. The EPO peptide of claim 10, wherein said fungal cell is a yeast cell.

12. A glycoPEGylated EPO peptide comprising an EPO peptide and at least one glycan and at least one poly(ethylene glycol) molecule covalently attached to said glycan,

wherein said poly(ethylene glycol) molecule is added to said EPO peptide using a glycosyltransferase.

13. The glycoPEGylated EPO peptide of claim 12, comprising at least one  
5 mono-antennary glycan.

14. The glycoPEGylated EPO peptide of claim 12, wherein all of said glycans  
are N-linked and are mono-antennary.

15. The glycoPEGylated EPO peptide of claim 12, wherein all of said glycans  
are N-linked and at least one of said glycans comprise said poly(ethylene glycol).

16. The glycoPEGylated EPO peptide of claim 15, wherein more than one of  
said glycans comprises said poly(ethylene glycol).

17. The glycoPEGylated EPO peptide of claim 12, wherein all of said glycans  
are N-linked and all of said glycans comprise said poly(ethylene glycol).

18. The glycoPEGylated EPO peptide of claim 12, comprising at least three  
20 mono-antennary glycans having said poly(ethylene glycol) covalently attached thereto.

19. A glycoPEGylated EPO peptide, wherein said EPO peptide comprises  
three or more glycans.

20. The glycoPEGylated EPO peptide of claim 9, wherein at least one of said  
glycans comprises said poly(ethylene glycol) covalently attached thereto.

21. The glycoPEGylated EPO peptide of claim 18, wherein more than one of  
said glycans comprises said poly(ethylene glycol) covalently attached thereto.

22. The glycoPEGylated EPO peptide of claim 18, wherein all of said glycans comprise said poly(ethylene glycol) covalently attached thereto.

23. The glycoPEGylated EPO peptide of claim 12 wherein said poly(ethylene glycol) is linked to at least one sugar moiety selected from the group consisting of fucose (Fuc), N-acetylglucosamine (GlcNAc), galactose (Gal) and a sialic acid (SA).

24. The glycoPEGylated EPO peptide of claim 23, wherein said sialic acid is N-acetylneuraminic acid.

25. The glycoPEGylated EPO peptide of claim 12, wherein said EPO peptide does not comprise an O-linked glycan.

26. The glycoPEGylated EPO peptide of claim 12 wherein said EPO peptide comprises at least one O-linked glycan.

27. The glycoPEGylated EPO peptide of claim 26, wherein said O-linked peptide comprises said poly(ethylene glycol) covalently attached thereto.

28. The glycoPEGylated EPO peptide of claim 27, wherein said EPO peptide is recombinantly expressed in a cell.

29. The glycoPEGylated EPO peptide of claim 28, wherein said cell is selected from the group consisting of an insect cell, a fungal cell and a mammalian cell.

30. The glycoPEGylated EPO peptide of claim 29, wherein said fungal cell is a yeast cell.

31. The glycoPEGylated EPO peptide of claim 29, wherein said cell is an insect cell.

32. The glycoPEGylated EPO peptide of claim 29, wherein said cell is a yeast cell.

33. The glycoPEGylated EPO peptide of claim 29, wherein said cell is a mammalian cell.

34. The glycoPEGylated EPO peptide of claim 33, wherein said mammalian cell is a CHO cell.

35. The glycoPEGylated EPO peptide of claim 12, wherein said poly(ethylene glycol) has a molecular weight selected from the group consisting of about 1 kDa, 2 kDa, 5 kDa, 10 kDa, 20 kDa, 30 kDa and 40 kDa.

36. The glycoPEGylated EPO peptide of claim 35, wherein said poly(ethylene glycol) has a molecular weight of 20 kDa.

37. The glycoPEGylated EPO peptide of claim 12, wherein said EPO peptide is selected from the group consisting of a naturally occurring EPO peptide and a mutated EPO peptide.

38. The glycoPEGylated EPO peptide of claim 37, wherein said mutated EPO peptide comprises the amino acid sequence of SEQ ID NO:73 having at least one mutation selected from the group consisting of Arg<sup>139</sup> to Ala<sup>139</sup>, Arg<sup>143</sup> to Ala<sup>143</sup> and Lys<sup>154</sup> to Ala<sup>154</sup>.

39. A method of making a glycoPEGylated EPO peptide, said method comprising the step of:

(a) contacting an EPO peptide with a mixture comprising a nucleotide sugar covalently linked to poly(ethylene glycol) and a glycosyltransferase under conditions sufficient to transfer said poly(ethylene glycol) to said EPO peptide.

40. The method of claim 39, wherein the sugar of said nucleotide sugar is selected from the group consisting of fucose (Fuc), N-acetylglucosamine (GlcNAc), galactose (Gal) and a sialic acid (SA).

5 41. The method of claim 40, wherein said sialic acid is N-acetylneuraminic acid (NAN).

42. The method of claim 39, wherein said poly(ethylene glycol) has a molecular weight selected from the group consisting of about 1 kDa, 2 kDa, 5 kDa, 10 kDa,  
10 20 kDa, 30 kDa and 40 kDa.

43. The method of claim 42, wherein said poly(ethylene glycol) has a molecular weight of 20 kDa.

15 44. The method of claim 39, wherein said EPO peptide is recombinantly expressed in a cell.

45. The method of claim 44, wherein said cell is selected from the group consisting of an insect cell, a fungal cell and a mammalian cell.

20 46. The method of claim 45, wherein said cell is an insect cell.

47. The method of claim 45, wherein said cell is a yeast cell.

25 48. The method of claim 45, wherein said cell is a mammalian cell.

49. The method of claim 48, wherein said mammalian cell is a CHO cell.

50. The method of claim 39, wherein said EPO peptide is selected from the  
30 group consisting of a naturally occurring EPO peptide and a mutated EPO peptide.

51. The method of claim 50, wherein said mature EPO peptide has the sequence of SEQ ID NO:73.

52. The method of claim 50, wherein said mutated EPO peptide comprises the  
5 amino acid sequence of SEQ ID NO: 73 having at least one mutation selected from the group consisting of Arg<sup>139</sup> to Ala<sup>139</sup>, Arg<sup>143</sup> to Ala<sup>143</sup> and Lys<sup>154</sup> to Ala<sup>154</sup>.

53. The method of claim 39, wherein before step (a):

(b) contacting said EPO peptide with a mixture comprising a nucleotide-N-  
10 acetylglucosamine (GlcNAc) molecule and an N-acetylglucosamine transferase (GnT) for which the nucleotide-GlcNAc is a substrate under conditions sufficient to form a bond between said GlcNAc and said EPO, wherein said GnT is selected from the group consisting of GnT I, GnT II, GnT III, GnT IV, GnT V and GnT VI.

54. The method of claim 53, wherein said mixture comprises one GnT  
15 selected from the group consisting of GnT I, GnT II, GnT IV, GnT V and GnT VI.

55. The method of claim 54, wherein said GnT is GnT I.

56. The method of claim 54, wherein said GnT is GnT II.

57. The method of claim 39, wherein said glycoPEGylated EPO peptide  
comprises at least one mono-antennary glycan.

58. The method of claim 39, wherein the sugar of said nucleotide sugar is  
25 galactose and said glycosyltransferase is galactosyl transferase I (GalT I).

59. The method of claim 53, wherein before step (a) but after step (b):

(c) contacting said EPO peptide with a mixture comprising a nucleotide galactose  
30 (Gal) and galactosyl transferase I (GalT I) under conditions sufficient to transfer galactose to said EPO peptide.

60. The method of claim 39, wherein in step (a), the sugar of said nucleotide sugar is sialic acid and said glycosyltransferase is a sialyltransferase.

5           61. The method of claim 60, wherein said sialic acid is N-acetylneuraminic acid (NAN).

62. The method of claim 60, wherein said sialyltransferase is selected from the group consisting of  $\alpha(2,3)$ sialyltransferase,  $\alpha(2,6)$ sialyltransferase and  
10   (2,8)sialyltransferase.

63. A glycoPEGylated EPO peptide made by the method of claim 39.

64. A glycoPEGylated EPO peptide, said EPO peptide comprising the  
15   sequence of SEQ ID NO:73.

65. A glycoPEGylated EPO peptide, said EPO peptide comprising the sequence of SEQ ID NO:73 and further comprising a mutation in said sequence.

20           66. A method of making a glycoPEGylated EPO peptide, said method comprising the steps of:

(a) contacting an EPO peptide with a mixture comprising a nucleotide sugar covalently linked to poly(ethylene glycol) and a glycosyltransferase under conditions sufficient to transfer said poly(ethylene glycol) to said EPO peptide, wherein said  
25   glycosyltransferase is a fucosyltransferase.

67. The method of claim 66, wherein said fucosyltransferase is selected from the group consisting of fucosyltransferase I, fucosyltransferase III, fucosyltransferase IV, fucosyltransferase V, fucosyltransferase VI and fucosyltransferase VII.

30           68. A glycoPEGylated EPO peptide made by the method of claim 66.

69. The method of claim 66, wherein said EPO peptide is expressed in a CHO cell.

70. A method of treating a mammal having anemia, said method comprising administering to said mammal an EPO peptide having one or more glycans having a glycoconjugate molecule attached to said peptide, wherein said EPO peptide is administered in an amount effective to increase the hematocrit level in said mammal.

71. The method of claim 70, wherein said mammal is a human.

72. A method of providing erythropoietin therapy to a mammal, said method comprising administering an effective amount of a glycoPEGylated EPO peptide comprising an EPO peptide and at least one glycan and at least one poly(ethylene glycol) molecule covalently attached to said glycan, wherein said poly(ethylene glycol) molecule is added to said EPO peptide using a glycosyltransferase, wherein said EPO peptide is administered in an amount effective to increase the hematocrit level in said mammal.

73. The method of claim 72, wherein said mammal is a human.

74. A method of treating a mammal having anemia, said method comprising administering to said mammal a glycoPEGylated EPO peptide comprising an EPO peptide and at least one glycan and at least one poly(ethylene glycol) molecule covalently attached to said glycan, wherein said poly(ethylene glycol) molecule is added to said EPO peptide using a glycosyltransferase, wherein said EPO peptide is administered in an amount effective to increase the hematocrit level in said mammal..

75. The method of claim 74, wherein said mammal is a human.



76. The method of claim 75, wherein said anemia is associated with chemotherapy.

77. A method of treating a kidney dialysis patient, said method comprising  
5 administering to said patient a glycoPEGylated EPO peptide comprising an EPO peptide and  
at least one glycan and at least one poly(ethylene glycol) molecule covalently attached to said  
glycan, wherein said poly(ethylene glycol) molecule is added to said EPO peptide using a  
glycosyltransferase, wherein said EPO peptide is administered in an amount effective to  
increase the hematocrit level in said patient.

10

1/498

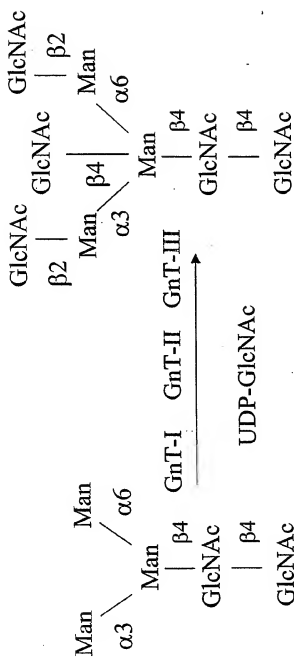
Trimannosyl core with  
Bisecting GlcNAc

FIG. 1

Trimannosyl core

2/498


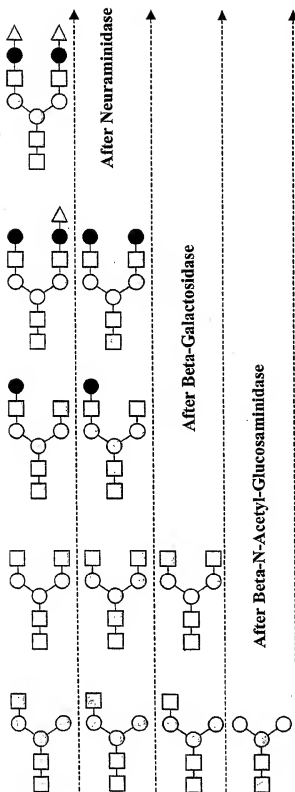


FIG. 2

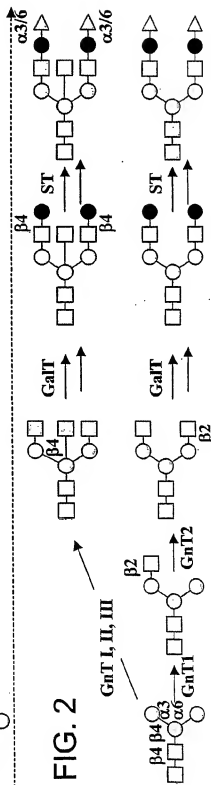

3/498

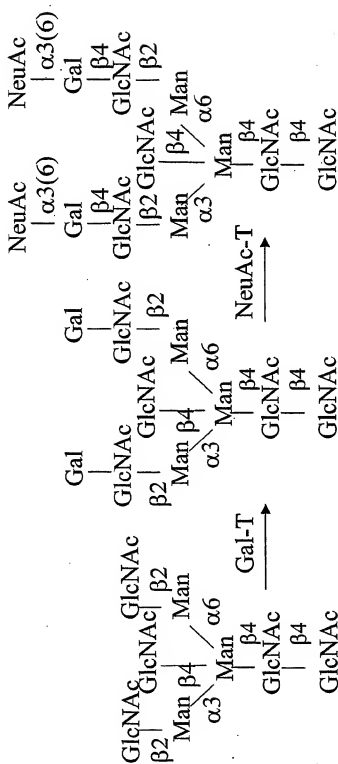


FIG. 3

4/498

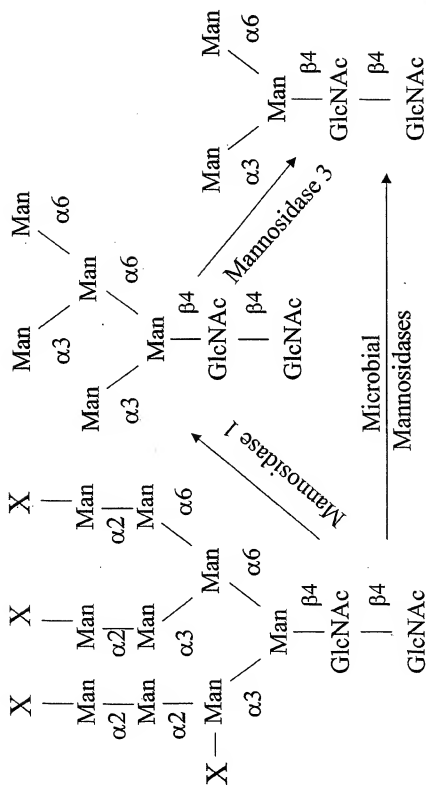


FIG. 4

5/498

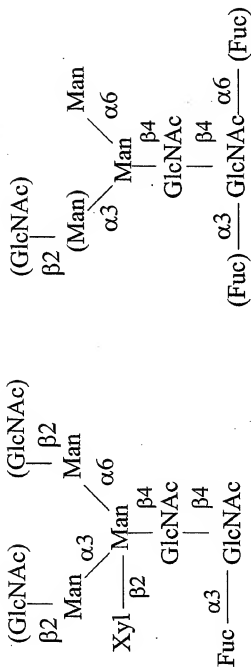


FIG. 6

FIG. 5



7/498

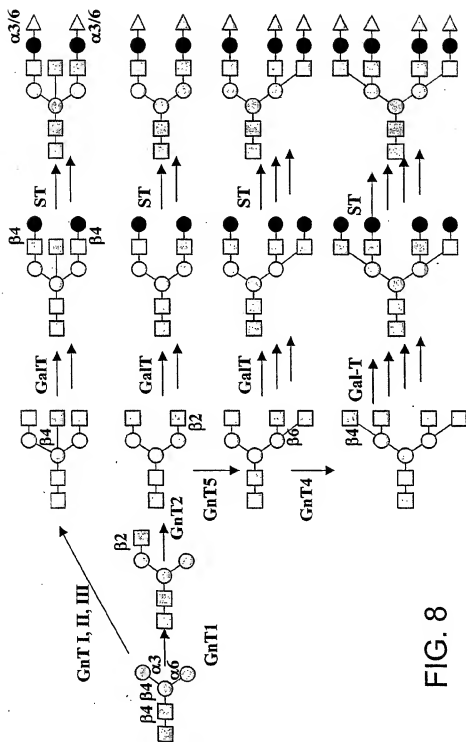


FIG. 8



8/498


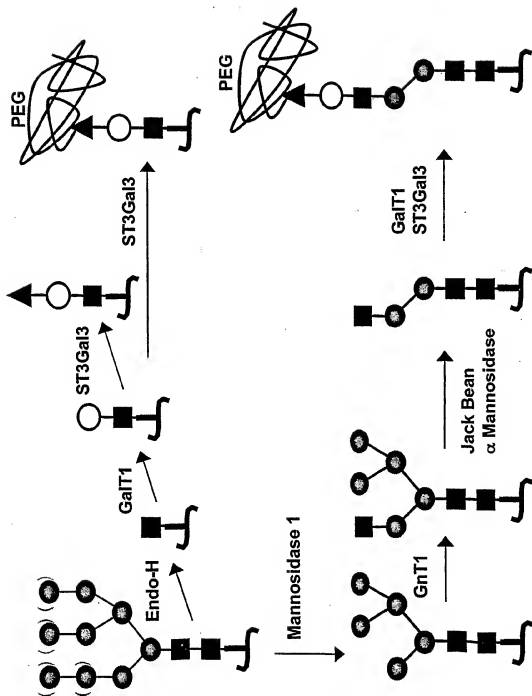


FIG. 9

9/498

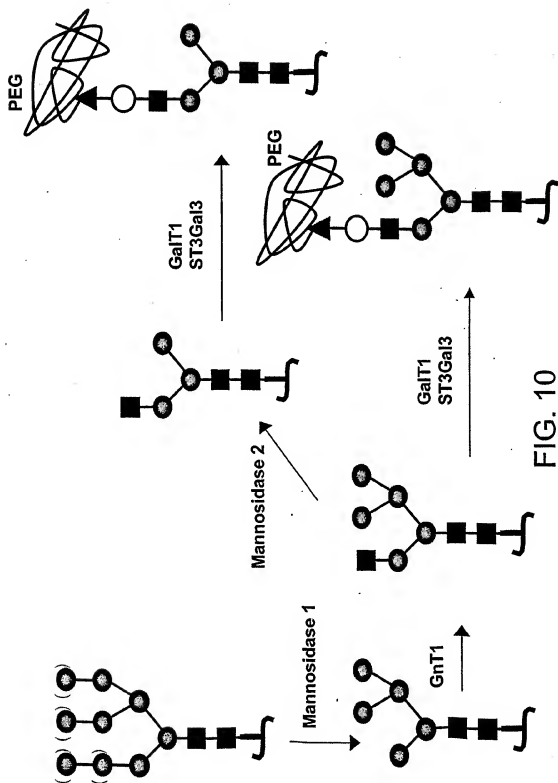

10/498

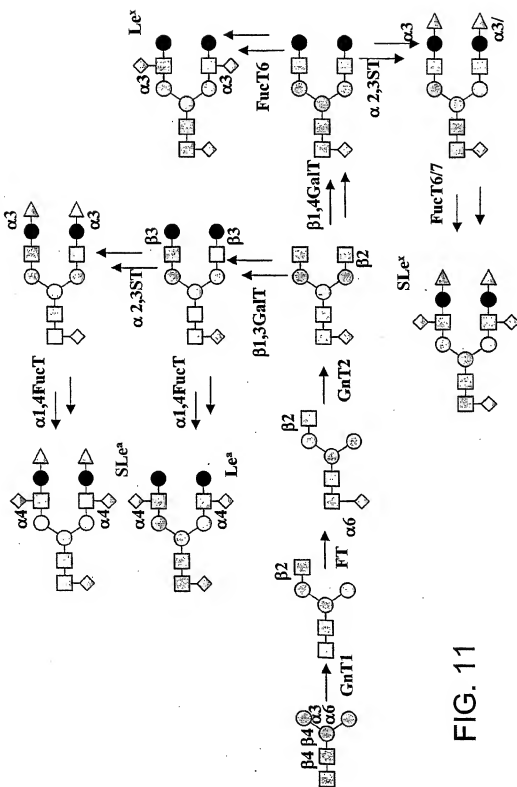


FIG. 11

11/498

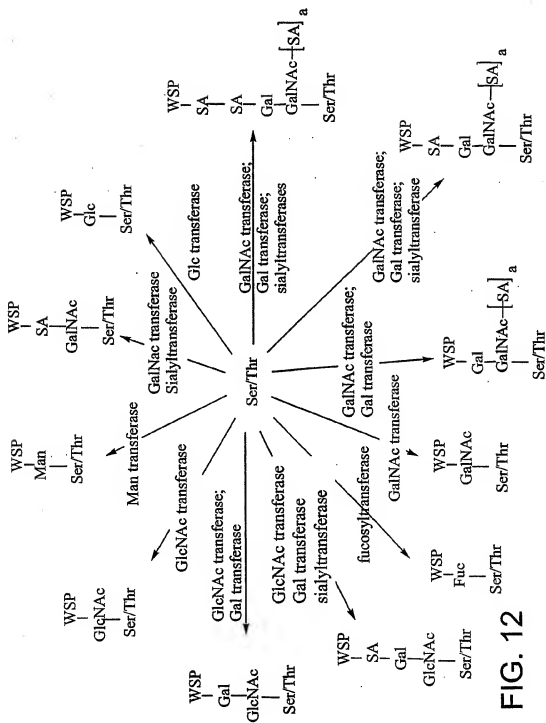


FIG. 12

12/498







**Core 1****Core 2****Core 3****Core 4**

FIG. 13



14/498



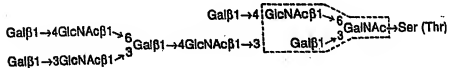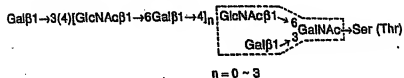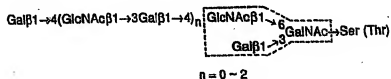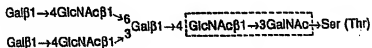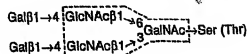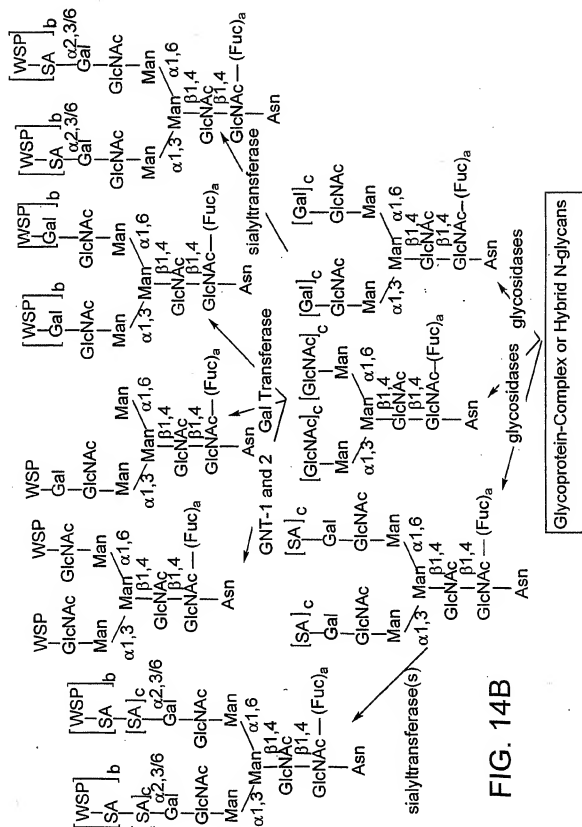















22/498


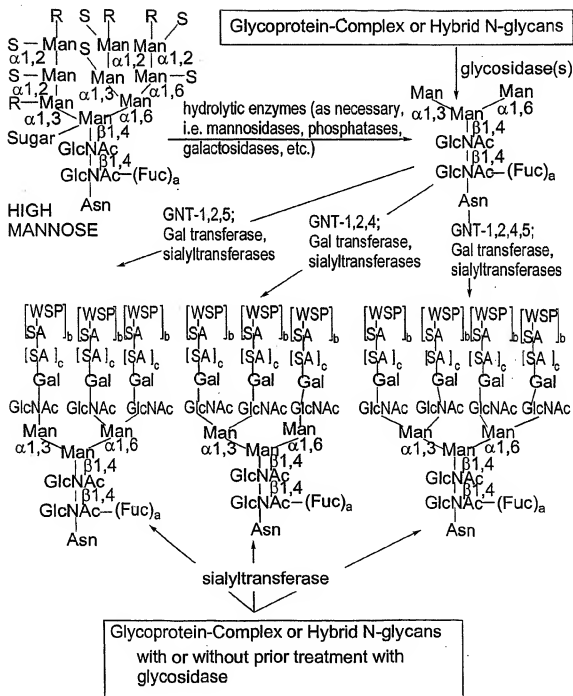


FIG. 20





24/498

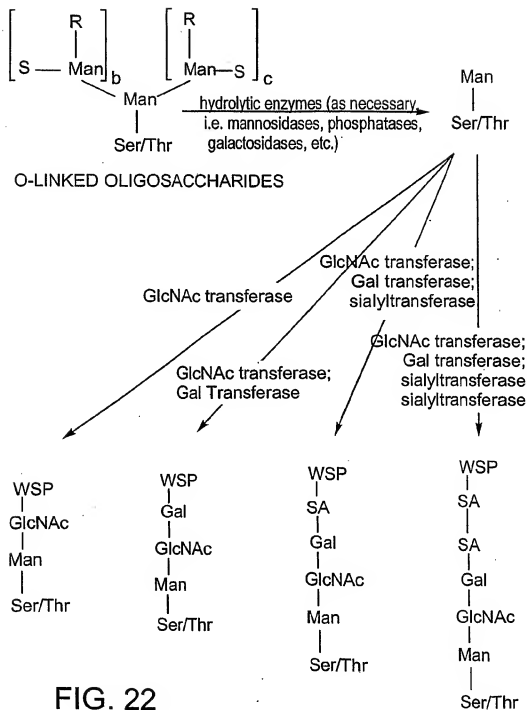

25/498

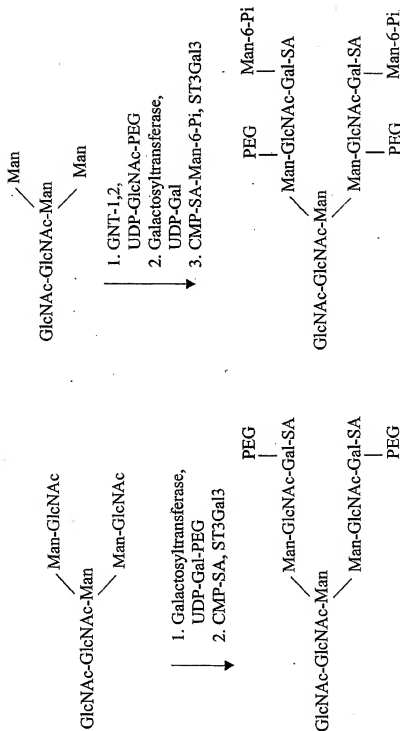


FIG. 23B

FIG. 23A

26/498

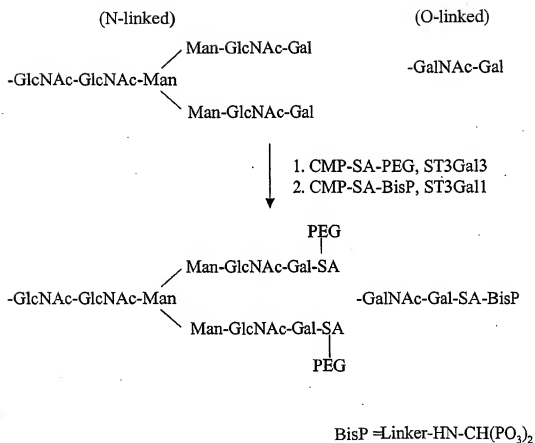


FIG. 23C

27/498

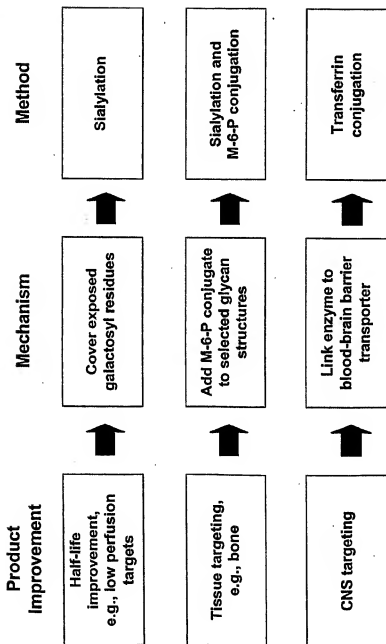


FIG. 24

28/498

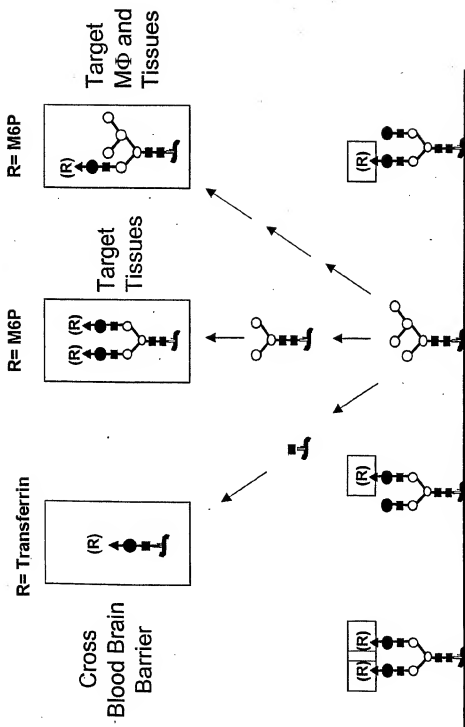


FIG. 25

29/498

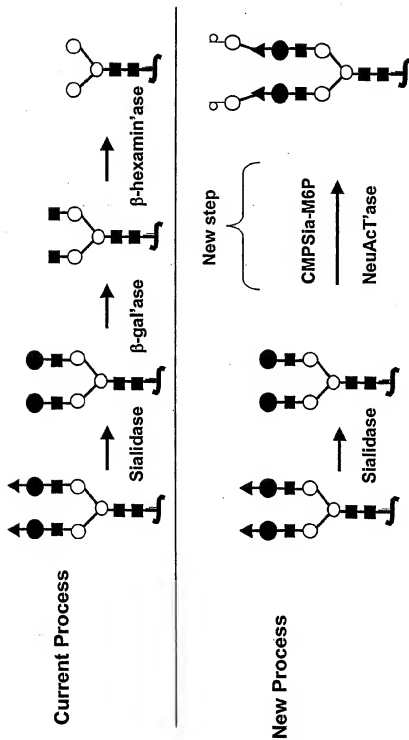


FIG. 26

30/498

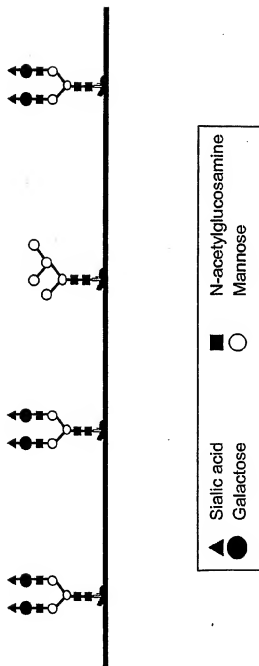


FIG. 27

## 31/498

12AP1/E5 -- Viventia Biotech  
 1964 -- Aventis  
 20K growth hormone -- AMUR  
 28P6/E6 -- Viventia Biotech  
 3-Hydroxyphthaloyl-beta-lactoglobulin --  
 4-IBB ligand gene therapy --  
 64-Cu MAb conjugate TETA-1A3 --  
 Mallinckrodt Institute of Radiology  
 64-Cu MAb conjugate TETA-cT84.66  
 64-Cu Trastuzumab TETA conjugate --  
 Genentech  
 A 200 -- Amgen  
 A10255 -- Eli Lilly  
 A1PDX -- Hedral Therapeutics  
 A6 -- Angstrom  
 aaAT-III -- Genzyme  
 Abciximab -- Centocor  
 ABI.001 -- Atlantic BioPharmaceuticals  
 ABT-828 -- Abbott  
 Accutin  
 Actinohivin  
 activin -- Biotech Australia, Human  
 Therapeutics, Curis  
 AD 439 -- Tanox  
 AD 519 -- Tanox  
 Adalimumab -- Cambridge Antibody Tech.  
 Adenocarcinoma vaccine -- Biomira -- NIS  
 Adenosine deaminase -- Enzond  
 Adenosine A2B receptor antagonists --  
 Adenosine Therapeutics  
 ADP-001 -- Axis Genetics  
 AF 13948 -- Affymax  
 Afelimomab -- Knoll  
 AFP-SCAN -- Immunomedics  
 AG 2195 -- Corixa  
 agalsidase alfa -- Transkaryotic Therapies  
 agalsidase beta -- Genzyme  
 AGENT-- Antisoma  
 AI 300 -- Autolimmune  
 AI-101 -- Teva  
 AI-102 -- Teva  
 AI-201 -- Autolimmune  
 AI-301 -- Autolimmune  
 AIDS vaccine -- ANRS, CIBG, Heseid  
 Biomed, Hollis-Eden, Rome, United  
 Biomedical, American Home Products,  
 Maxygen  
 airway receptor ligand -- IC Innovations  
 AJW 2 -- Ajinomoto  
 AK 30 NGF -- Alkermes  
 Albuferon -- Human Genome Sciences  
 albumin -- Biogen, DSM Anti-Infectives,  
 Genzyme Transgenics, PPL Therapeutics,  
 TranXenoGen, Welfide Corp.  
 aldesleukin -- Chiron  
 alefacept -- Biogen  
 Alemtuzumab  
 Allergy therapy -- ALK-Abello/Maxygen,  
 ALK-Abello/RP Scherer  
 allergy vaccines -- Allergy Therapeutics  
 Alnidofibatide -- Aventis Pasteur  
 Alnorine -- SRC VB VECTOR  
 ALP 242 -- Gruenenthal  
 Alpha antitrypsin -- Arriva/Hyland  
 Immuno/ProMetic/Protease Sciences  
 Alpha-1 antitrypsin -- Cutter, Bayer, PPL  
 Therapeutics, Profile, ZymoGenetics,  
 Arriva  
 Alpha-1 protease inhibitor -- Genzyme  
 Transgenics, Welfide Corp.  
 Alpha-galactose fusion protein --  
 Immunomedics  
 Alpha-galactosidase A -- Research  
 Corporation Technologies, Genzyme  
 Alpha-glucosidase -- Genzyme, Novazyme  
 Alpha-lactalbumin  
 Alpha-L-iduronidase -- Transkaryotic  
 Therapies, BioMarin  
 alteplase -- Genentech  
 alvircept sudotox -- NIH  
 ALX-0600, a GLP-2 agonist -- NPS Allelix  
 Corp.

FIG. 28A



32/498

ALX1-11 --sNPS Pharmaceuticals  
 Alzheimer's disease gene therapy  
 AM-133 -- AMRAD  
 Amb a 1 immunostim conj. -- Dynavax  
 AMD 3100 -- AnorMED -- NIS  
 AMD 3465 -- AnorMED -- NIS  
 AMD 3465 -- AnorMED -- NIS  
 AMD Fab -- Genentech  
 Amediplase -- Menarini, Novartis  
 AM-F9  
 Amoebiasis vaccine  
 Amphiregulin -- Octagene  
 anakinra -- Amgen  
 analgesic -- Nobex  
 aneastim -- Amgen  
 Anergix.RA -- Corixa, Organon  
 Angiocidin -- InKine  
 angiogenesis inhibitors -- ILEX  
 AngioMab -- Antisoma  
 Angiopoietins -- Regeneron/Procter &  
 Gamble  
 angiostatin -- EntreMed  
 Angiostatin/endostatin gene therapy --  
 Genetix Pharmaceuticals  
 angiotensin-II, topical -- Maret  
 Anthrax -- EluSys Therapeutics/US Army  
 Medical Research Institute  
 Anthrax vaccine  
 Anti platelet-derived growth factor D human  
 monoclonal antibodies -- CuraGen  
 Anti-17-1A Mab 3622W94 --  
 GlaxoSmithKline  
 Anti-2C4 Mab -- Genentech  
 anti-4-1BB monoclonal antibodies -- Bristol-  
 Myers Squibb  
 Anti-Adhesion Platform Tech. -- CytoVax  
 Anti-adipocyte Mab -- Cambridge Antibody  
 Tech./ObeSys  
 antiallergics -- Maxygen  
 antiallergy vaccine -- Acambis  
 Anti-alpha-4-integrin Mab  
 Anti-alphavβ3 integrin Mab -- Applied  
 Molecular Evolution  
 Anti-angiogenesis monoclonal antibodies --  
 KS Biomedix/Schering AG  
 Anti-B4 Mab-DC1 conjugate -- ImmunoGen  
 Anti-B7 antibody PRIMATIZED -- IDEC  
 Anti-B7-1 Mab 16-10A1  
 Anti-B7-1 Mab 1G10  
 Anti-B7-2 Mab GL-1  
 Anti-B7-2-gelonin immunotoxin --  
 Antibacterials/antifungals --  
 Diversa/IntraBiotics  
 Anti-beta-amyloid monoclonal antibodies --  
 Cambridge Antibody Tech., Wyeth-Ayerst  
 Anti-BLyS antibodies -- Cambridge  
 Antibody Tech./Human Genome Sciences  
 Antibody-drug conjugates -- Seattle  
 Genetics/Eos  
 Anti-C5 Mab BB5-1 -- Alexion  
 Anti-C5 Mab N19-8 -- Alexion  
 Anti-C8 Mab  
 anticancer cytokines -- BioPulse  
 anticancer matrix -- Telios Integra  
 Anticancer monoclonal antibodies -- ARIUS,  
 Immunex  
 anticancer peptides -- Maxygen, Micrologix  
 Anticancer prodrug Tech. -- Alexion  
 Antibody Technologies  
 anticancer Troy-Bodies -- Affite -- Affitech  
 anticancer vaccine -- NIH  
 anticancers -- Epimmune  
 Anti-CCR5/CXCR4 sheep Mab -- KS  
 Biomedix Holdings  
 Anti-CD11a Mab KBA --  
 Anti-CD11a Mab M17  
 Anti-CD11a Mab TA-3 --  
 Anti-CD11a Mab WT.1 --  
 Anti-CD11b Mab -- Pharmacia  
 Anti-CD11b Mab LM2  
 Anti-CD154 Mab -- Biogen  
 Anti-CD16-anti-CD30 Mab -- Biotest

FIG. 28B

## 33/498

Anti-CD18 MAb -- Pharmacia	Anti-CD4 MAb -- Centocor, IDEC
Anti-CD19 MAb B43 --	Pharmaceuticals, Xenova Group
Anti-CD19 MAb -liposomal sodium butyrate conjugate --	Anti-CD4 MAb 16H5
Anti-CD147	Anti-CD4 MAb 4162W94 -- GlaxoSmithKline
Anti-CD19 MAb-saporin conjugate --	Anti-CD4 MAb B-F5 -- Diaclone
Anti-CD19-dsFv-PE38-immunotoxin --	Anti-CD4 MAb GK1-5
Anti-CD2 MAb 12-15 --	Anti-CD4 MAb KT6
Anti-CD2 MAb B-E2 -- Diaclone	Anti-CD4 MAb OX38
Anti-CD2 MAb OX34 --	Anti-CD4 MAb PAP conjugate -- Bristol-Myers Squibb
Anti-CD2 MAb OX54 --	Anti-CD4 MAb RIB 5-2
Anti-CD2 MAb OX55 --	Anti-CD4 MAb W3/25
Anti-CD2 MAb RM2-1	Anti-CD4 MAb YTA 3.1.2
Anti-CD2 MAb RM2-2	Anti-CD4 MAb YTS 177-9
Anti-CD2 MAb RM2-4	Anti-CD40 ligand MAb 5c8 -- Biogen
Anti-CD20 MAb BCA B20	Anti-CD40 MAb
Anti-CD20-anti-Fc alpha RI bispecific MAb -- Medarex, Tenovus	Anti-CD40 MAb 5D12 -- Tanox
Anti-CD22 MAb-saporin-6 complex --	Anti-CD44 MAb A3D8
Anti-CD3 Immunotoxin --	Anti-CD44 MAb GKWA3
Anti-CD3 MAb 145-2C11 -- Pharming	Anti-CD44 MAb IM7
Anti-CD3 MAb CD4IgG conjugate -- Genentech	Anti-CD44 MAb KM81
Anti-CD3 MAb humanised -- Protein Design, RW Johnson	Anti-CD44 variant monoclonal antibodies -- Corixa/Hebrew University
Anti-CD3 MAb WT32	Anti-CD45 MAb BC8-I-131
Anti-CD3 MAb-ricin-chain-A conjugate --	Anti-CD45RB MAb
Anti-CD3 MAb-xanthine-oxidase conjugate --	Anti-CD48 MAb HuLy-m3
Anti-CD30 MAb BerH2 -- Medac	Anti-CD48 MAb WM-63
Anti-CD30 MAb-saporin conjugate	Anti-CD5 MAb -- Becton Dickinson
Anti-CD30-scFv-ETA'-immunotoxin	Anti-CD5 MAb OX19
Anti-CD38 MAb AT13/5	Anti-CD6 MAb
Anti-CD38 MAb-saporin conjugate	Anti-CD7 MAb-PAP conjugate
Anti-CD3-anti-CD19 bispecific MAb	Anti-CD7 MAb-ricin-chain-A conjugate
Anti-CD3-anti-EGFR MAb	Anti-CD8 MAb -- Amerimmune, Cytodyn, Becton Dickinson
Anti-CD3-anti-interleukin-2-receptor MAb	Anti-CD8 MAb 2-43
Anti-CD3-anti-MOV18 MAb -- Centocor	Anti-CD8 MAb OX8
Anti-CD3-anti-SCLC bispecific MAb	Anti-CD80 MAb P16C10 -- IDEC
Anti-CD4 idiotype vaccine	Anti-CD80 MAb P7C10 -- ID Vaccine
	Anti-CD8-idarubicin conjugate
	Anti-CEA MAb CE-25
	Anti-CEA MAb MN 14 -- Immunomedics

FIG. 28C

## 34/498

- Anti-CEA MAb MN14-PE40 conjugate -- Immunomedics
- Anti-CEA MAb T84.66-interleukin-2 conjugate
- Anti-CEA sheep MAb -- KS Biomedix Holdings
- Anti-cell surface monoclonal antibodies -- Cambridge Antibody Tech. /Pharmacia
- Anti-c-erbB2-anti-CD3 bifunctional MAb -- Otsuka
- Anti-CMV MAb -- Scotgen
- Anti-complement
- Anti-CTLA-4 MAb
- Anti-EGFR catalytic antibody -- Hersed Biomed
- anti-EGFR immunotoxin -- IVAX
- Anti-EGFR MAb -- Abgenix
- Anti-EGFR MAb 528
- Anti-EGFR MAb KSB 107 -- KS Biomedix
- Anti-EGFR MAb-DM1 conjugate -- ImmunoGen
- Anti-EGFR MAb-LA1 --
- Anti-EGFR sheep MAb -- KS Biomedix
- Anti-FAP MAb F19-I-131
- Anti-Fas IgM MAb CH11
- Anti-Fas MAb Jo2
- Anti-Fas MAb RK-8
- Anti-Fit-1 monoclonal antibodies -- ImClone
- Anti-fungal peptides -- State University of New York
- antifungal tripeptides -- BTG
- Anti-ganglioside GD2 antibody-interleukin-2 fusion protein -- Lexigen
- Anti-GM2 MAb -- Kyowa
- Anti-GM-CSF receptor monoclonal antibodies -- AMRAD
- Anti-gp130 MAb -- Tosoh
- Anti-HCA monoclonal antibodies -- AltaRex/Epigen
- Anti-hCG antibodies -- Abgenix/AVI BioPharma
- Anti-heparanase human monoclonal antibodies -- Oxford Glycosciences/Medarex
- Anti-hepatitis C virus human monoclonal antibodies -- XTL Biopharmaceuticals
- Anti-HER-2 antibody gene therapy
- Anti-herpes antibody -- Epicyte
- Anti-HIV antibody -- Epicyte
- anti-HIV catalytic antibody -- Hersed Biomed
- anti-HIV fusion protein -- Idun
- anti-HIV proteins -- Cangene
- Anti-HM1-24 MAb -- Chugai
- Anti-hrR3 MAb
- Anti-Human-Carcinoma-Antigen MAb -- Epicyte
- Anti-ICAM-1 MAb -- Boehringer Ingelheim
- Anti-ICAM-1 MAb 1A-29 -- Pharmacia
- Anti-ICAM-1 MAb HA58
- Anti-ICAM-1 MAb YN1/1.7.4
- Anti-ICAM-3 MAb ICM3 -- ICOS
- Anti-idiotypic breast cancer vaccine 11D10
- Anti-idiotypic breast cancer vaccine ACA14C5 --
- Anti-idiotypic cancer vaccine -- ImClone Systems/Merck KGaA ImClone, Viventia Biotech
- Anti-idiotypic cancer vaccine 1A7 -- Titan
- Anti-idiotypic cancer vaccine 3H1 -- Titan
- Anti-idiotypic cancer vaccine TriAb -- Titan
- Anti-idiotypic Chlamydia trachomatis vaccine
- Anti-idiotypic colorectal cancer vaccine -- Novartis
- Anti-idiotypic colorectal cancer vaccine -- Onyvox
- Anti-idiotypic melanoma vaccine -- IDEC Pharmaceuticals
- Anti-idiotypic ovarian cancer vaccine ACA 125
- Anti-idiotypic ovarian cancer vaccine AR54 - AltaRex

FIG. 28D

## 35/498

Anti-idiotype ovarian cancer vaccine CA-125 – AltaRex, Biomira	Anti-L-selectin monoclonal antibodies -- Protein Design Labs, Abgenix, Stanford University
Anti-IgE catalytic antibody -- Hesed Biomed	Anti-MBL monoclonal antibodies -- Alexion/Brigham and Women's Hospital
Anti-IgE MAb E26 -- Genentech	Anti-MHC monoclonal antibodies
Anti-IGF-1 MAb	Anti-MIF antibody humanised -- IDEC, Cytokine PharmaSciences
anti-inflammatory -- GeneMax	Anti-MRSA/VRSA sheep MAb -- KS Biomedix Holdings
anti-inflammatory peptide -- BTG	Anti-mu MAb -- Novartis
anti-integrin peptides -- Burnha	Anti-MUC-1 MAb
Anti-interferon-alpha-receptor MAb 64G12 -- Pharma Pacific Management	Anti-MUC 18
Anti-interferon-gamma MAb -- Protein Design Labs	Anti-Nogo-A MAb IN1
Anti-interferon-gamma polyclonal antibody - Advanced Biotherapy	Anti-nuclear autoantibodies -- Procyon
Anti-interleukin-10 MAb --	Anti-ovarian cancer monoclonal antibodies -
Anti-interleukin-12 MAb --	- Dompe
Anti-interleukin-1-beta polyclonal antibody -- R&D Systems	Anti-p185 monoclonal antibodies
Anti-interleukin-2 receptor MAb 2A3	Anti-p43 MAb
Anti-interleukin-2 receptor MAb 33B3-1 -- Immunotech	Antiparasitic vaccines
Anti-interleukin-2 receptor MAb ART-18	Anti-PDGF/bFGF sheep MAb -- KS Biomedix
Anti-interleukin-2 receptor MAb LO-Tact-1	Anti-properdin monoclonal antibodies -- Abgenix/Gliatech
Anti-interleukin-2 receptor MAb Mikbeta1	Anti-PSMA (prostrate specific membrane antigen)
Anti-interleukin-2 receptor MAb NDS61	Anti-PSMA MAb J591 -- BZL Biologics
Anti-interleukin-4 MAb 11B11	Anti-Rev MAb gene therapy --
Anti-interleukin-5 MAb -- Wallace Laboratories	Anti-RSV antibodies -- Epicyte, Intracell
Anti-interleukin-6 MAb -- Centocor, Diacorne, Pharmadigm	Anti-RSV monoclonal antibodies -- Medarex/MedImmune, Applied Molecular Evolution/MedImmune
Anti-interleukin-8 MAb -- Abgenix	Anti-RSV MAb, inhalation -- Alkermes/MedImmune
Anti-interleukin-8 MAb -- Xenotech	Anti-RT gene therapy
Anti-JL1 MAb	Antisense K-ras RNA gene therapy
Anti-Klebsiella sheep MAb -- KS Biomedix Holdings	Anti-SF-25 MAb
Anti-Laminin receptor MAb-liposomal doxorubicin conjugate	Anti-sperm antibody -- Epicyte
Anti-LCG MAb -- Cytoclonal	Anti-Tac(Fv)-PE38 conjugate
Anti-lipopolysaccharide MAb -- VitaResc	Anti-TAPA/CD81 MAb AMP1
	Anti-tat gene therapy

36/498

Anti-TCR-alphabeta MAb H57-597  
 Anti-TCR-alphabeta MAb R73  
 Anti-tenascin MAb BC-4-I-131  
 Anti-TGF-beta human monoclonal  
 antibodies -- Cambridge Antibody Tech.,  
 Genzyme  
 Anti-TGF-beta MAb 2G7 -- Genentech  
 Antithrombin III -- Genzyme Transgenics,  
 Aventis, Bayer, Behringwerke, CSL,  
 Myriad  
 Anti-Thy1 MAb  
 Anti-Thy1.1 MAb  
 Anti-tissue factor/factor VIIA sheep MAb --  
 KS Biomedix  
 Anti-TNF monoclonal antibodies --  
 Centocor, Chiron, Peptech, Pharacia,  
 Sero  
 Anti-TNF sheep MAb -- KS Biomedix  
 Holdings  
 Anti-TNFalpha MAb -- Genzyme  
 Anti-TNFalpha MAb B-C7 -- Diaclone  
 Anti-tooth decay MAb -- Planet BioTech.  
 Anti-TRAIL receptor-1 MAb -- Takeda  
 Antitumour RNases -- NIH  
 Anti-VCAM MAb 2A2 -- Alexion  
 Anti-VCAM MAb 3F4 -- Alexion  
 Anti-VCAM-1 MAb  
 Anti-VEC MAb -- ImClone  
 Anti-VEGF MAb -- Genentech  
 Anti-VEGF MAb 2C3  
 Anti-VEGF sheep MAb -- KS Biomedix  
 Holdings  
 Anti-VLA-4 MAb HP1/2 -- Biogen  
 Anti-VLA-4 MAb PS/2  
 Anti-VLA-4 MAb R1-2  
 Anti-VLA-4 MAb TA-2  
 Anti-VAP-1 human MAb  
 Anti-VRE sheep MAb -- KS Biomedix  
 Holdings  
 ANUP -- TranXenoGen  
 ANUP-1 -- Pharis

AOP-RANTES -- Senetek  
 Apan-CH -- Praecis Pharmaceuticals  
 APC-8024 -- Demegen  
 ApoA-1 -- Milano, Pharmacia  
 Apogen -- Alexion  
 apolipoprotein A1 -- Avanir  
 Apolipoprotein E -- Bio-Tech. General  
 Applaggin -- Biogen  
 aprotinin -- ProdiGene  
 APT-070C -- AdProTech  
 AR 177 -- Aronex Pharmaceuticals  
 AR 209 -- Aronex Pharmaceuticals,  
 Antigenics  
 AR545C  
 ARGENT gene delivery systems -- ARIAD  
 Arresten  
 ART-123 -- Asahi Kasei  
 arylsulfatase B -- BioMarin  
 arylsulfatase B, Recombinant human --  
 BioMarin  
 AS 1051 -- Ajinomoto  
 ASI-BCL -- Intracell  
 Asparaginase - Merck  
 ATL-101 -- Alizyme  
 Atrial natriuretic peptide -- Pharis  
 Aurintricarboxylic acid-high molecular  
 weight  
 Autoimmune disorders -- GPC  
 Biotech/MorphoSys  
 Autoimmune disorders and transplant  
 rejection -- Bristol-Myers Squibb/Genzyme  
 Tra  
 Autoimmune disorders/cancer --  
 Abgenix/Chiron, CuraGen  
 Autotaxin  
 Avicidin -- NeoRx  
 axogenesis factor-1 -- Boston Life Sciences  
 Axokine -- Regeneron  
 B cell lymphoma vaccine -- Biomira  
 B7-1 gene therapy --  
 BABS proteins -- Chiron

FIG. 28F

## 37/498

BAM-002 -- Novelos Therapeutics  
 Basiliximab (anti CD25 MAb) -- Novartis  
 Bay-16-9996 -- Bayer  
 Bay-39-9437 -- Bayer  
 Bay-50-4798 -- Bayer  
 BB-10153 -- British Biotech  
 BBT-001 -- Bolder BioTech.  
 BBT-002 -- Bolder BioTech.  
 BBT-003 -- Bolder BioTech.  
 BBT-004 -- Bolder BioTech.  
 BBT-005 -- Bolder BioTech.  
 BBT-006 -- Bolder BioTech.  
 BBT-007 -- Bolder BioTech.  
 BCH-2763 -- Shire  
 BCSF -- Millenium Biologix  
 BDNF -- Regeneron -- Amgen  
 Becapiermin -- Johnson & Johnson, Chiron  
 Bectumomab -- Immunomedics  
 Beriplast -- Aventis  
 Beta-adrenergic receptor gene therapy --  
 University of Arkansas  
 bFGF -- Scios  
 BI 51013 -- Behringwerke AG  
 BIBH 1 -- Boehringer Ingelheim  
 BIM-23190 -- Beaufour-Ipsen  
 birch pollen immunotherapy -- Pharmacia  
 bispecific fusion proteins -- NIH  
 Bispecific MAb 2B1 -- Chiron  
 Bitistatin  
 BIWA 4 -- Boehringer Ingelheim  
 blood substitute -- Northfield, Baxter Intl.  
 BLP-25 -- Biomira  
 BLS-0597 -- Boston Life Sciences  
 BLyS -- Human Genome Sciences  
 BLyS radiolabelled -- Human Genome  
 Sciences  
 BM 06021 -- Boehringer Mannheim  
 BM-202 -- BioMarin  
 BM-301 -- BioMarin  
 BM-301 -- BioMarin  
 BM-302 -- BioMarin  
 BMP 2 -- Genetics Institute/Medtronic-  
 Sofamor Danek, Genetics Institute/  
 Collagenesis, Genetics  
 Institute/Yamanouch  
 BMP 2 gene therapy  
 BMP 52 -- Aventis Pasteur, Biopharm  
 BMP-2 -- Genetics Institute  
 BMS 182248 -- Bristol-Myers Squibb  
 BMS 202448 -- Bristol-Myers Squibb  
 bone growth factors -- IsoTis  
 BPC-15 -- Pfizer  
 brain natriuretic peptide --  
 Breast cancer -- Oxford  
 GlycoSciences/Medarex  
 Breast cancer vaccine -- Therion Biologics,  
 Oregon  
 BSSL -- PPL Therapeutics  
 BST-2001 -- BioStratum  
 BST-3002 -- BioStratum  
 BTI 322 --  
 butyrylcholinesterase -- Shire  
 C 6822 -- COR Therapeutics  
 C1 esterase inhibitor -- Pharming  
 C3d adjuvant -- AdProTech  
 CAB-2.1 -- Millennium  
 calcitonin -- Inhale Therapeutics Systems,  
 Aventis, Genetronics, TranXenoGen,  
 Unigene, Rhone Poulenc Rohrer  
 calcitonin -- oral -- Nobex, Emisphere,  
 Pharmaceutical Discovery  
 Calcitonin gene-related peptide -- Asahi  
 Kasei -- Unigene  
 calcitonin, human -- Suntory  
 calcitonin, nasal -- Novartis, Unigene  
 calcitonin, Panoderm -- Elan  
 calcitonin, Peptitrol -- Shire  
 calcitonin, salmon -- Therapicon  
 calin -- Biopharm  
 Calphobindin I  
 calphobindin I -- Kowa  
 calreticulin -- NYU

FIG. 28G

## 38/498

Campath-1G  
 Campath-1M  
 cancer therapy -- Cangene  
 cancer vaccine -- Aixlie, Aventis Pasteur,  
 Center of Molecular Immunology, YM  
 BioSciences, Cytos, Genzyme,  
 Transgenics, GlobelImmune, Igeneon,  
 ImClone, Virogenetics, InterCell, Iomai,  
 Jenner Biotherapies, Memorial Sloan-  
 Kettering Cancer Center, Sydney Kimmel  
 Cancer Center, Novavax, Protein  
 Sciences, Argonex, SIGA  
 Cancer vaccine ALVAC-CEA B7.1 --  
 Aventis Pasteur/Therion Biologics  
 Cancer vaccine CEA-TRICOM -- Aventis  
 Pasteur/Therion Biologics  
 Cancer vaccine gene therapy -- Cantab  
 Pharmaceuticals  
 Cancer vaccine HER-2/neu -- Corixa  
 Cancer vaccine THERATOPE -- Biomira  
 cancer vaccine, PolyMASC -- Valentis  
 Candida vaccine -- Corixa, Inhibitex  
 Canstatin -- ILEX  
 CAP-18 -- Panorama  
 Cardiovascular gene therapy -- Collateral  
 Therapeutics  
 carperitide -- Suntory  
 Casocidin-1 -- Pharis  
 CAT 152 -- Cambridge Antibody Tech.  
 CAT 192 -- Cambridge Antibody Tech.  
 CAT 213 -- Cambridge Antibody Tech.  
 Catalase-- Enzon  
 Cat-PAD -- Circassia  
 CB 0006 -- Celltech  
 CCK(27-32)-- Akzo Nobel  
 CCR2-64I -- NIH  
 CD, Procept -- Paligent  
 CD154 gene therapy  
 CD39 -- Immunex  
 CD39-L2 -- Hyseq  
 CD39-L4 -- Hyseq  
 CD4 fusion toxin -- Senetek  
 CD4 IgG -- Genentech  
 CD4 receptor antagonists --  
 Pharmacocepeia/Progenics  
 CD4 soluble -- Progenics  
 CD4, soluble -- Genzyme Transgenics  
 CD40 ligand -- Immunex  
 CD4-ricin chain A -- Genentech  
 CD59 gene therapy -- Alexion  
 CD8 TIL cell therapy -- Aventis Pasteur  
 CD8, soluble -- Avidex  
 CD95 ligand -- Roche  
 CDP 571 -- Celltech  
 CDP 850 -- Celltech  
 CDP-860 (PEG-PDGF MAb) -- Celltech  
 CDP 870 -- Celltech  
 CDS-1 -- Ernest Orlando  
 Cedelizumab -- Ortho-McNeil  
 Cetermin -- Insmad  
 CETP vaccine -- Avant  
 Cetorelix  
 Cetuximab  
 CGH 400 -- Novartis  
 CGP 42934 -- Novartis  
 CGP 51901 -- Tanox  
 CGRP -- Unigene  
 CGS 27913 -- Novartis  
 CGS 32359 -- Novartis  
 Chagas disease vaccine -- Corixa  
 chemokines -- Immune Response  
 CHH 380 -- Novartis  
 chitinase -- Genzyme, ICOS  
 Chlamydia pneumoniae vaccine -- Antex  
 Biologics  
 Chlamydia trachomatis vaccine -- Antex  
 Biologics  
 Chlamydia vaccine -- GlaxoSmithKline  
 Cholera vaccine CVD 103-HgR -- Swiss  
 Serum and Vaccine Institute Berne  
 Cholera vaccine CVD 112 -- Swiss Serum  
 and Vaccine Institute Berne

FIG. 28H

## 39/498

Cholera vaccine inactivated oral -- SBL	CRL 1605 -- CytRx
Vaccin	CS-560 -- Sankyo
Chrysalin -- Chrysalis BioTech.	CSF -- ZymoGenetics
CI-782 -- Hitachi Kase	CSF-G -- Hangzhou, Dong-A, Hanmi
Ciliary neurotrophic factor -- Fidia, Roche	CSF-GM -- Cangene, Hunan, LG Chem
CIM project -- Active Biotech	CSF-M -- Zarix
CL 329753 -- Wyeth-Ayerst	CT 1579 -- Merck Frosst
CL22, Cobra -- ML Laboratories	CT 1786 -- Merck Frosst
Clenoliximab -- IDEC	CT-112 <sup>A</sup> -- BTG
Clostridium difficile antibodies -- Epicycle	CTB-134L -- Xenova
clotting factors -- Octagene	CTC-111 -- Kaketsuken
CMB 401 -- Celltech	CTGF -- FibroGen
CNTF -- Sigma-Tau	CTLA4-Ig -- Bristol-Myers Squibb
Cocaine abuse vaccine -- Cantab,	CTLA4-Ig gene therapy --
ImmuLogic, Scripps	CTP-37 -- AVI BioPharma
coccidiomycosis vaccine -- Arizo	C-type natriuretic peptide -- Suntory
collagen -- Type I -- Pharming	CVS 995 -- Corvas Intl.
Collagen formation inhibitors -- FibroGen	CX 397 -- Nikko Kyodo
Collagen/hydroxyapatite/bone growth factor	CY 1747 -- Epimmune
-- Aventis Pasteur, Biopharm, Orquest	CY 1748 -- Epimmune
collagenase -- BioSpecifics	Cyanovirin-N
Colorectal cancer vaccine -- Wistar Institute	Cystic fibrosis therapy -- CBR/IVAX
Component B, Recombinant -- Sero	CYT 351
Connective tissue growth factor inhibitors --	cytokine Traps -- Regeneron
FibroGen/Taisho	cytokines -- Enzon, Cytodonal
Contortostat	Cytomegalovirus glycoprotein vaccine --
contraceptive vaccine -- Zonagen	Chiron, Aquila Biopharmaceuticals,
Contraceptive vaccine hCG	Aventis Pasteur, Virogenetics
Contraceptive vaccine male reversible --	Cytomegalovirus vaccine live -- Aventis
IMMUCON	Pasteur
Contraceptive vaccine zona pellucida --	Cytosine deaminase gene therapy --
Zonagen	GlaxoSmithKline
Copper-64 labelled Mab TETA-1A3 -- NCI	DA-3003 -- Dong-A
Coralyne	DAB389interleukin-6 -- Senetek
Corsevin M	DAB389interleukin-7
C-peptide analogues -- Schwarz	DAC:GLP-2 -- ConjuChem, Inc.
CPI-1500 -- Consensus	Daclizumab (anti-IL2R Mab) -- Protein
CRF -- Neurobiological Tech.	Design Labs
cRGDfV pentapeptide --	DAMP <sup>A</sup> -- Incyte Genomics
CRL 1095 -- CytRx	Daniplestim -- Pharmacia
CRL 1336 -- CytRx	darbepoetin alfa -- Amgen



40/498

DBI-3019 -- Diabetogen  
 DCC -- Genzyme  
 DDF -- Hyseq  
 decorin -- Integra, Telios  
 defensins -- Large Scale Biology  
 DEGR-VIIa  
 Delimmunised antibody 3B6/22 AGEN  
 Deimmunised anti-cancer antibodies --  
     Biovation/Viragen  
 Dendroamide A  
 Dengue vaccine -- Bavarian Nordic, Merck  
 denileukin diftotox -- Ligand  
 DES-1101 -- Desmos  
 desirudin -- Novartis  
 desmopressin -- Unigene  
 Desmoteplase -- Merck, Schering AG  
 Destablase  
 Diabetes gene therapy -- DeveloGen, Pfizer  
 Diabetes therapy -- Crucell  
 Diabetes type 1 vaccine -- Diamyd  
     Therapeutics  
 DiaCIM -- YM BioSciences  
 dialytic oligopeptides -- Research Corp  
 Diamyd -- Diamyd Therapeutics  
 DiaPep227 -- Peppen  
 DiavaX -- Corixa  
 Digoxin MAb -- Glaxo  
 Diphtheria tetanus pertussis-hepatitis B  
     vaccine -- GlaxoSmithKline  
 DIR therapy -- Solis Therapeutics --  
 DNase -- Genentech  
 Dornase alfa -- Genentech  
 Dornase alfa, inhalation -- Genentech  
 Doxorubicin-anti-CEA MAb conjugate --  
     Immunomedics  
 DP-107 -- Trimeris  
 drotrecogin alfa -- Eli Lilly  
 DTctGMCSF  
 DTP-polio vaccine -- Aventis Pasteur  
 DU 257-KM231 antibody conjugate --  
     Kyowa  
     dural graft matrix -- Integra  
     Duteplase -- Baxter Intl.  
     DWP-401 -- Daewoong  
     DWP-404 -- Daewoong  
     DWP-408 -- Daewoong  
     Dx 88 (Epi-KAL2) -- Dyax  
     Dx 890 (elastin inhibitors) -- Dyax  
     E coli O157 vaccine -- NIH  
     E21-R -- BresaGen  
     Eastern equine encephalitis virus vaccine --  
     Echicetin --  
     Echinhibin 1 --  
     Echistatin -- Merck  
     Echitamine --  
     Ecromeximab -- Kyowa Hakko  
     EC-SOD -- PPL Therapeutics  
     Eculizumab (5G1.1) -- Alexion  
     EDF -- Ajinomoto  
     EDN derivative -- NIH  
     EDNA -- NIH  
     Edobacomab -- XOMA  
     Edrecolomab -- Centocor  
     EF 5077  
     Efalizumab -- Genentech  
     EGF fusion toxin -- Seragen, Ligand  
     EGF-P64k vaccine -- Center of Molecular  
         Immunology  
     EL 246 -- LigoCyte  
     elastase inhibitor -- Synergen  
     elcatonin -- Therapicon  
     EMD 72000 -- Merck KGaA  
     Emdogain -- BIORA  
     emfilermin -- AMRAD  
     Emoctakin -- Novartis  
     enamel matrix protein -- BIORA  
     Endo III -- NYU  
     endostatin -- EntreMed, Pharis  
     Enhancins -- Micrologix  
     Enlimomab -- Isis Pharm.  
     Enoxaparin sodium -- Phamuka

FIG. 28J

41/498

enzyme linked antibody nutrient depletion  
 therapy -- KS Biomedix Holdings  
 Eosinophil-derived neutralizing agent --  
 EP-51216 -- Asta Medica  
 EP-51389 -- Asta Medica  
 EPH family ligands -- Regeneron  
 Epidermal growth factor -- Hitachi Kasei,  
 Johnson & Johnson  
 Epidermal growth factor fusion toxin --  
 Senetek  
 Epidermal growth factor-genistein --  
 EPI-HNE-4 -- Dyax  
 EPI-KAL2 -- Dyax  
 Epoetin-alfa -- Amgen, Dragon  
 Pharmaceuticals, Nanjing Huaxin  
 Epratuzumab -- Immunomedics  
 Epstein-Barr virus vaccine --  
 Aviron/SmithKline Beecham, Bioresearch  
 Eptacog alfa -- Novo Nordisk  
 Eptifibatide -- COR Therapeutics  
 erb-38 --  
 Erlizumab -- Genentech  
 erythropoietin -- Alkermes, ProLease, Dong-  
 A, Elanex, Genetics Institute, LG Chem,  
 Protein Sciences, Serono, Snow Brand,  
 SRC VB VECTOR, Transkaryotic  
 Therapies  
 Erythropoietin Beta -- Hoffman La Roche  
 Erythropoietin/Epoetin alfa -- Chugai  
 Escherichia coli vaccine -- North American  
 Vaccine, SBL Vaccin, Swiss Serum and  
 Vaccine Institute Berne  
 etanercept -- Immunex  
 examorelin -- Mediolanum  
 Exendin 4 -- Amylin  
 exonuclease VII  
 F 105 -- Centocor  
 F-992 -- Fornix  
 Factor IX -- Alpha Therapeutics, Welfide  
 Corp., CSL, enetics Institute/AHP,  
 Pharmacia, PPL Therapeutics  
 Factor IX gene therapy -- Cell Genesys  
 Factor VII -- Novo Nordisk, Bayer, Baxter  
 Intl.  
 Factor VIIa -- PPL Therapeutics,  
 ZymoGenetics  
 Factor VIII -- Bayer Genentech, Beaufour-  
 Ipsen, CLB, Inex, Octagen, Pharmacia,  
 Pharming  
 Factor VIII -- PEGylated -- Bayer  
 Factor VIII fragments -- Pharmacia  
 Factor VIII gene therapy -- Targeted  
 Genetics  
 Factor VIII sucrose formulation -- Bayer,  
 Genentech  
 Factor VIII-2 -- Bayer  
 Factor VIII-3 -- Bayer  
 Factor Xa inhibitors -- Merck, Novo Nordisk,  
 Mochida  
 Factor XIII -- ZymoGenetics  
 Factors VIII and IX gene therapy -- Genetics  
 Institute/Targeted Genetics  
 Famoxin -- Genset  
 Fas (delta) TM protein -- LXR BioTech.  
 Fas TR -- Human Genome Sciences  
 Felvizumab -- Scotgen  
 FFR-VIIa -- Novo Nordisk  
 FG-001 -- F-Gene  
 FG-002 -- F-Gene  
 FG-004 -- F-Gene  
 FG-005 -- F-Gene  
 FGF + fibrin -- Repair  
 Fibrimage -- Bio-Tech. General  
 fibrin-binding peptides -- ISIS Innovation  
 fibrinogen -- PPL Therapeutics, Pharming  
 fibroblast growth factor -- Chiron, NYU,  
 Ramot, ZymoGenetics  
 fibrolase conjugate -- Schering AG  
 Filgrastim -- Amgen  
 filgrastim -- PDA modified -- Xencor  
 FLT-3 ligand -- Immunex  
 FN18 CRM9 --

FIG. 28K

## 42/498

follostatin -- Biotech Australia, Human Therapeutics	Glucocerebrosidase -- Genzyme
folliotropin alfa -- Alkermes, ProLease, PowderJect, Serono, Akzo Nobel	glutamate decarboxylase -- Genzyme Transgenics
Follitropin Beta -- Bayer, Organon FP 59	Glycoprotein S3 -- Kureha
FSH -- Ferring	GM-CSF -- Immunhex
FSH + LH -- Ferring	GM-CSF tumour vaccine -- PowderJect
F-spondin -- CeNeS	GnRH immunotherapeutic -- Protherics
fusion protein delivery system -- UAB Research Foundation	Goserelin (LhRH antagonist) -- AstraZeneca
fusion toxins -- Boston Life Sciences	gp75 antigen -- ImClone
G 5598 -- Genentech	gp96 -- Antigenics
GA-II -- Transkaryotic Therapies	GPI 0100 -- Galenica
Gamma-interferon analogues -- SRC VB VECTOR	GR 4991W93 -- GlaxoSmithKline
Ganirelix -- Roche	Granulocyte colony-stimulating factor -- Dong-A
gastric lipase -- Meristem	Granulocyte colony-stimulating factor conjugate
Gavilimomab --	grass allergy therapy -- Dynavax
G-CSF -- Amgen, SRC VB VECTOR	GRF1-44 -- ICN
GDF-1 -- CeNeS	Growth Factor -- Chiron, Atrigel, Atrix, Innogenetics, ZymoGenetics, Novo
GDF-5 -- Biopharm	growth factor peptides -- Biotherapeutics
GDNF (glial derived neurotrophic factor) -- Amgen	growth hormone -- LG Chem
gelsolin -- Biogen	growth hormone, Recombinant human -- Serono
Gemtuzumab ozogamicin -- Celltech	GT 4086 -- Gliatech
Gene-activated epoetin-alfa -- Aventis Pharma -- Transkaryotic Therapies	GW 353430 -- GlaxoSmithKline
Glanzmann thrombasthenia gene therapy --	GW-278884 -- GlaxoSmithKline
Glatiramer acetate -- Yeda	H 11 -- Viventia Biotech
glial growth factor 2 -- CeNeS	H5N1 influenza A virus vaccine -- Protein Sciences
GLP-1 -- Amylin, Suntory, TheraTech, Watson	haemoglobin -- Biopure
GLP-1 peptide analogues -- Zealand Pharmaceuticals	haemoglobin 3011, Recombinant -- Baxter Healthcare
GLP-2 -- Novo Nordisk, Ontario, Inc., Suntory Limited	haemoglobin crosumaril -- Baxter Intl.
glucagon -- Eli Lilly, ZymoGenetics	haemoglobin stabilized -- Ajinomoto
Glucagon-like peptide-1 7-36 amide -- Suntory	haemoglobin, recombinant -- Apex
Glucogen-like peptide -- Amylin	HAF -- Immune Response
	Hantavirus vaccine
	HB 19
	HBNF -- Regeneron
	HCC-1 -- Pharis

FIG. 28L

## 43/498

hCG -- Milkhaus  
 hCG vaccine -- Zonagen  
 HE-317 -- Hollis-Eden Pharmaceuticals  
 Heat shock protein cancer and influenza vaccines -- StressGen  
 Helicobacter pylori vaccine -- Acambis, AstraZeneca/CSL, Chiron, Provalis  
 Helistat-G -- GalaGen  
 Hemolink -- Hemosol  
 hepapoietin -- Snow Brand  
 heparanase -- InSight  
 heparinase I -- Ibex  
 heparinase III -- Ibex  
 Hepatitis A vaccine -- American Biogenetic Sciences  
 Hepatitis A vaccine inactivated  
 Hepatitis A vaccine Nothav -- Chiron  
 Hepatitis A-hepatitis B vaccine -- GlaxoSmithKline  
 hepatitis B therapy -- Tripep  
 Hepatitis B vaccine -- Amgen, Chiron SpA, Meiji Milk, NIS, Prodeva, PowderJect, Rhein Biotech  
 Hepatitis B vaccine recombinant -- Evans Vaccines, Epitex Combiotech, Genentech, MedImmune, Merck Sharp & Dohme, Rhein Biotech, Shantha Biotechnics, Vector, Yeda  
 Hepatitis B vaccine recombinant TGP 943 -- Takeda  
 Hepatitis C vaccine -- Bavarian Nordic, Chiron, Innogenetics Acambis,  
 Hepatitis D vaccine -- Chiron Vaccines  
 Hepatitis E vaccine recombinant -- Genelabs/GlaxoSmithKline, Novavax  
 hepatocyte growth factor -- Panorama, Sosei  
 hepatocyte growth factor kringle fragments -  
 - Entremed  
 Her-2/Neu peptides -- Corixa  
 Herpes simplex glycoprotein DNA vaccine -- Merck, Wyeth-Lederle Vaccines-Malvern, Genentech, GlaxoSmithKline, Chiron, Takeda  
 Herpes simplex vaccine -- Cantab Pharmaceuticals, CEL-SCI, Henderson Morley  
 Herpes simplex vaccine live -- ImClone Systems/Wyeth-Lederle, Aventis Pasteur  
 HGF derivatives -- Dompe  
 hAPP vaccine -- Crucell  
 Hib-hepatitis B vaccine -- Aventis Pasteur  
 HIC 1  
 HIP -- Altachem  
 Hirudins -- Biopharma, Cangene, Dongkook, Japan Energy Corporation, Pharmacia Corporation, SIR International, Sanofi-Synthelabo, Sotragene, Rhein Biotech  
 HIV edible vaccine -- ProdiGene  
 HIV gp120 vaccine -- Chiron, Ajinomoto, GlaxoSmithKline, ID Vaccine, Progenics, VaxGen  
 HIV gp120 vaccine gene therapy --  
 HIV gp160 DNA vaccine -- PowderJect, Aventis Pasteur, Oncogen, Hyland Immuno, Protein Sciences  
 HIV gp41 vaccine -- Panacos  
 HIV HGP-30W vaccine -- CEL-SCI  
 HIV immune globulin -- Abbott, Chiron  
 HIV peptides -- American Home Products  
 HIV vaccine -- Applied bioTech., Axis Genetics, Biogen, Bristol-Myers Squibb, Genentech, Korea Green Cross, NIS, Oncogen, Protein Sciences Corporation, Terumo, Tonen Corporation, Wyeth-Ayerst, Wyeth-Lederle Vaccines-Malvern, Advanced BioScience Laboratories, Bavarian Nordic, Bavarian Nordic/Statens Serum Institute, GeneCure, Immune Response, Progenics, Theron Biologics, United Biomedical, Chiron

FIG. 28M

## 44/498

HIV vaccine vCP1433 -- Aventis Pasteur  
 HIV vaccine vCP1452 -- Aventis Pasteur  
 HIV vaccine vCP205 -- Aventis Pasteur  
 HL-9 -- American BioScience  
 HM-9239 -- Cytran  
 HML-103 -- Hemosol  
 HML-104 -- Hemosol  
 HML-105 -- Hemosol  
 HML-109 -- Hemosol  
 HML-110 -- Hemosol  
 HML-121 -- Hemosol  
 hNLP -- Pharis  
 Hookworm vaccine  
 host-vector vaccines -- Henogen  
 HPM 1 -- Chugai  
 HPV vaccine -- MediGene  
 HSA -- Meristem  
 HSF -- StressGen  
 HSP carriers -- Weizmann, Yeda, Peptor  
 HSPPC-70 -- Antigenics  
 HSPPC-96, pathogen-derived -- Antigenics  
 HSV 863 -- Novartis  
 HTLV-I DNA vaccine  
 HTLV-I vaccine  
 HTLV-II vaccine -- Access  
 HU 901 -- Tanox  
 Hu23F2G -- ICOS  
 HuHMF1  
 HumaLYM -- Intracell  
 Human krebs statika -- Yamanouchi  
 human monoclonal antibodies --  
   Abgenix/Biogen, Abgenix/ Corixa,  
   Abgenix/immunex, Abgenix/Lexicon,  
   Abgenix/ Pfizer, Athersys/Medarex,  
   Biogen/MorphoSys, CAT/Searle,  
   Centocor/Medarex, Corixa/Kirin Brewery,  
   Corixa/Medarex, Eos BioTech./Medarex,  
   Eos/Xenex, Exelixis/Protein Design  
   Labs, ImmunoGen/ Raven, Medarex/  
   B.Twelve, MorphoSys/ImmunoGen, XTL  
   Biopharmaceuticals/Dyax,  
 Human monoclonal antibodies --  
   Medarex/Northwest Biotherapeutics,  
   Medarex/Seattle Genetics  
 human netrin-1 -- Exelixis  
 human papillomavirus antibodies -- Epicyte  
 Human papillomavirus vaccine -- Biotech  
   Australia, IDEC, StressGen  
 Human papillomavirus vaccine MEDI 501 --  
   MedImmune/GlaxoSmithKline  
 Human papillomavirus vaccine MEDI  
   503/MEDI 504 --  
   MedImmune/GlaxoSmithKline  
 Human papillomavirus vaccine TA-CIN --  
   Cantab Pharmaceuticals  
 Human papillomavirus vaccine TA-HPV --  
   Cantab Pharmaceuticals  
 Human papillomavirus vaccine TH-GW --  
   Cantab/GlaxoSmithKline  
 human polyclonal antibodies -- Biosite/Eos  
   BioTech./ Medarex  
 human type II anti factor VIII monoclonal  
   antibodies -- ThromboGenics  
 humanised anti glycoprotein Ib murine  
   monoclonal antibodies -- ThromboGenics  
 HumaRAD -- Intracell  
 HuMax EGFR -- Genmab  
 HuMax-CD4 -- Medarex  
 HuMax-IL15 -- Genmab  
 HYB 190 -- Hybridon  
 HYB 676 -- Hybridon  
 I-125 Mab A33 -- Celltech  
 Ibritumomab tiuxetan -- IDEC  
 IBT-9401 -- Ibx  
 IBT-9402 -- Ibx  
 IC 14 -- ICOS  
 Idarubicin anti-Ly-2.1 --  
 IDEC 114 -- IDEC  
 IDEC 131 -- IDEC  
 IDEC 152 -- IDEC  
 IDM 1 -- IDM  
 IDPS -- Hollis-Eden Pharmaceuticals

FIG. 28N

## 45/498

iduronate-2-sulfatase -- Transkaryotic Therapies  
 IGF/IBP-2-13 -- Pharis  
 IGN-101 -- Igeneon  
 IK HIR02 -- Iketon  
 IL-11 -- Genetics Institute/AHP  
 IL-13-PE38 -- NeoPharm  
 IL-17 receptor -- Immunex  
 IL-18BP -- Yeda  
 IL-1Hy1 -- Hyseq  
 IL-1 $\beta$  -- Celltech  
 IL-1 $\beta$  adjuvant -- Celltech  
 IL-2 -- Chiron  
 IL-2 + IL-12 -- Hoffman La-Roche  
 IL-6/sIL-6R fusion -- Hadasit  
 IL-6R derivative -- Tosoh  
 IL-7-Dap 389 fusion toxin -- Ligand  
 IL-21 -- Novo Nordisk, ZymoGenetics  
 IM-862 -- Cytran  
 IMC-1C11 -- ImClone  
 imiglucerase -- Genzyme  
 Immune globulin intravenous (human) -- Hoffman La Roche  
 immune privilege factor -- Proneuron  
 Immunocal -- Immunotec  
 Immunogene therapy -- Briana Bio-Tech  
 Immunoliposomal 5-fluorodeoxyuridine-dipalmitate --  
 immunosuppressant vaccine -- Aixlie  
 immunotoxin -- Antisoma, NIH  
 ImmuRAIT-Re-188 -- Immunomedics  
 imreg-1 -- Imreg  
 infertility -- Johnson & Johnson, E-TRANS  
 Infliximab -- Centocor  
 Influenza virus vaccine -- Aventis Pasteur, Protein Sciences  
 inhibin -- Biotech Australia, Human Therapeutics  
 Inhibitory G protein gene therapy  
 INKP-2001 -- InKine  
 Inolimomab -- Diaclone  
 insulin -- AutoImmune, Altea, Biobras, BioSante, Bio-Tech. General, Chong Kun Dang, Emisphere, Flamel, Provalis, Rhein Biotech, TranXenoGen  
 insulin (bovine) -- Novartis  
 insulin analogue -- Eli Lilly  
 Insulin Aspart -- Novo Nordisk  
 insulin detemir -- Novo Nordisk  
 insulin glargine -- Aventis  
 insulin inhaled -- Inhale Therapeutics Systems, Alkermes  
 insulin oral -- Inovax  
 insulin, AeroDose -- AeroGen  
 insulin, AERx -- Aradigm  
 insulin, BEODAS -- Elan  
 insulin, Biphasix -- Helix  
 insulin, buccal -- Generex  
 insulin, I2R -- Flemington  
 insulin, intranasal -- Bentley  
 insulin, oral -- Nobex, Unigene  
 insulin, Orasome -- Endorex  
 insulin, ProMaxx -- Epic  
 insulin, Quadrant -- Elan  
 insulin, recombinant -- Aventis  
 insulin, Spiros -- Elan  
 insulin, Transfersome -- IDEA  
 insulin, Zymo, recombinant -- Novo Nordisk  
 insulinotropin -- Scios  
 Insulysin gene therapy --  
 integrin antagonists -- Merck  
 interferon (Alpha2) -- SRC VB VECTOR, Viragen, Dong-A, Hoffman La-Roche, Genentech  
 interferon -- BioMedicines, Human Genome Sciences  
 interferon (Alfa-n3) -- Interferon Sciences Intl.  
 interferon (Alpha), Biphasix -- Helix

FIG. 280

## 46/498

interferon (Alpha)—Amgen, BioNative,  
 Novartis, Genzyme Transgenics,  
 Hayashibara, Inhale Therapeutics  
 Systems, Medusa, Flamel, Dong-A,  
 GeneTrol, Nasteck, Shantha,  
 Wassermann, LG Chem, Sumitomo,  
 Aventis, Behring EGIS, Pepgen, Servier,  
 Rhein Biotech,  
 interferon (Alpha2A)  
 interferon (Alpha2B) — Enzon, Schering-  
 Plough, Biogen, IDEA  
 interferon (Alpha-N1) — GlaxoSmithKline  
 interferon (beta) — Rentschler, GeneTrol,  
 Meristem, Rhein Biotech, Toray, Yeda,  
 Daiichi, Mochida  
 interferon (Beta1A) — Sero, Biogen  
 interferon (beta1A), inhale — Biogen  
 interferon (β1b)— Chiron  
 interferon (tau)— Pepgen  
 Interferon alfacon-1 — Amgen  
 Interferon alpha-2a vaccine  
 Interferon Beta 1b — Schering/Chiron,  
 InterMune  
 Interferon Gamma — Boehringer Ingelheim,  
 Sheffield, Rentschler, Hayashibara  
 interferon receptor, Type I — Sero  
 interferon (Gamma1B) — Genentech  
 Interferon-alpha-2b + ribavirin — Biogen,  
 ICN  
 Interferon-alpha-2b gene therapy —  
 Schering-Plough  
 Interferon-con1 gene therapy —  
 interleukin-1 antagonists — Dompe  
 Interleukin-1 receptor antagonist — Abbott  
 Bioresearch, Pharmacia  
 Interleukin-1 receptor type I — Immunex  
 interleukin-1 receptor Type II — Immunex  
 Interleukin-1 trap — Regeneron  
 Interleukin-1-alpha — Immunex/Roche  
 interleukin-2 — SRC VB VECTOR,  
 Ajinomoto, Biomira, Chiron  
 IL-2/ diphtheria toxin — Ligand  
 Interleukin-3 — Cangene  
 Interleukin-4 — Immunology Ventures,  
 Sanofi Winthrop, Schering-Plough,  
 Immunex/ Sanofi Winthrop, Bayer, Ono  
 interleukin-4 + TNF-Alpha — NIH  
 interleukin-4 agonist — Bayer  
 interleukin-4 fusion toxin — Ligand  
 Interleukin-4 receptor — Immunex, Immun  
 Interleukin-6 — Ajinomoto, Cangene, Yeda,  
 Genetics Institute, Novartis  
 interleukin-6 fusion protein  
 interleukin-6 fusion toxin — Ligand, Sero  
 interleukin-7 — IC Innovations  
 interleukin-7 receptor — Immunex  
 interleukin-8 antagonists — Kyowa  
 Hakko/Millennium/Pfizer  
 interleukin-9 antagonists — Genaera  
 Interleukin-10 — DNAX, Schering-Plough  
 Interleukin-10 gene therapy —  
 interleukin-12 — Genetics Institute, Hoffman  
 La-Roche  
 interleukin-13 — Sanofi  
 interleukin-13 antagonists — AMRAD  
 Interleukin-13-PE38QQR  
 interleukin-15 — Immunex  
 interleukin-16 — Research Corp  
 Interleukin-18 — GlaxoSmithKline  
 Interleukin-18 binding protein — Sero  
 Ior-P3 — Center of Molecular Immunology  
 IP-10 — NIH  
 IPF — Metabolex  
 IR-501 — Immune Response  
 ISIS 9125 — Isis Pharmaceuticals  
 ISURF No. 1554 — Millennium  
 ISURF No. 1866 — Iowa State Univer.  
 ITF-1697 — Italfarmaco  
 IxC 162 — Ixion  
 J 695 — Cambridge Antibody Tech.,  
 Genetics Inst., Knoll  
 Jagged + FGF — Repair

## FIG. 28P

47/498

JKC-362 -- Phoenix Pharmaceuticals  
 JTP-2942 -- Japan Tobacco  
 Juman monoclonal antibodies --  
     Medarex/Raven  
 K02 -- Axyx Pharmaceuticals  
 Keliximab -- IDEC  
 Keyhole limpet haemocyanin  
 KGF -- Amgen  
 KM 871 -- Kyowa  
 KPI 135 -- Scios  
 KPI-022 -- Scios  
 Kringle 5  
 KSB 304  
 KSB-201 -- KS Biomedex  
 L 696418 -- Merck  
 L 703801 -- Merck  
 L1 -- Acorda  
 L-761191 -- Merck  
 lactoferrin -- Meristem, Pharming, Agennix  
 lactoferrin cardio -- Pharming  
 LAG-3 -- Seroxo  
 LAIT -- GEMMA  
 LAK cell cytotoxin -- Arizona  
 lamellarins -- PharmaMar/University of  
     Malaga  
 laminin A peptides -- NIH  
 lanoteplase -- Genetics Institute  
 laronidase -- BioMarin  
 Lassa fever vaccine  
 LCAT -- NIH  
 LDP 01 -- Millennium  
 LDP 02 -- Millennium  
 Lecithinized superoxide dismutase --  
     Seikagaku  
 LeIF adjuvant -- Corixa  
 leishmaniasis vaccine -- Corixa  
 lenercept -- Hoffman La-Roche  
 Lenograstim -- Aventis, Chugai  
 lepirudin -- Aventis  
 leptin -- Amgen, IC Innovations  
 Leptin gene therapy -- Chiron Corporation  
 leptin, 2nd-generation -- Amgen  
 leridistim -- Pharmacia  
 leuprolide, ProMaxx -- Epic  
 leuprorelin, oral -- Unigene  
 LeuTech -- Papatin  
 LEX 032 -- SuperGen  
 LiDEPT -- Novartis  
 Lintuzumab (anti-CD33 MAb) -- Protein  
     Design Labs  
 lipase -- Altus Biologics  
 lipid A vaccine -- EntreMed  
 lipid-linked anchor Tech. -- ICRT, ID  
     Biomedical  
 liposome-CD4 Tech. -- Sheffield  
 Listeria monocytogenes vaccine  
 LMB 1  
 LMB 7  
 LMB 9 -- Battelle Memorial Institute, NIH  
 LM-CD45 -- Cantab Pharmaceuticals  
 lovastatin -- Merck  
 LSA-3  
 LT- $\beta$  receptor -- Biogen  
 lung cancer vaccine -- Corixa  
 lusupultide -- Scios  
 L-Vax -- AVAX  
 LY 355455 -- Eli Lilly  
 LY 366405 -- Eli Lilly  
 LY-355101 -- Eli Lilly  
 Lyme disease DNA vaccine -- Vical/Aventis  
     Pasteur  
 Lyme disease vaccine -- Aquila  
 Biopharmaceuticals, Aventis, Pasteur,  
 Symbicom, GlaxoSmithKline, Hyland  
 Immuno, MedImmune  
 Lymphocytic choriomeningitis virus vaccine  
 lymphoma vaccine -- Biomira, Genitope  
 LYP18  
 lys plasminogen, recombinant  
 Lysosomal storage disease gene therapy --  
     Avigen  
 lysostaphin -- Nutrition 21

FIG. 28Q



## 48/498

M 23 -- Gruenenthal  
 M1 monoclonal antibodies -- Acorda  
 Therapeutics  
 MA 16N7C2 -- Corvas Intl.  
 malaria vaccine -- GlaxoSmithKline,  
 AdProTech, Antigenics, Apovia, Aventis  
 Pasteur, Axis Genetics, Behringwerke,  
 CDCP, Chiron Vaccines, Genzyme  
 Transgenics, Hawaii, MedImmune, NIH,  
 NYU, Oxxon, Roche/Saramane, Biotech  
 Australia, Rx Tech  
 Malaria vaccine CDC/NIIMALVAC-1  
 malaria vaccine, multicomponent  
 mammaglobin -- Corixa  
 mammastatin -- Biotherapeutics  
 mannan-binding lectin -- NatlImmu  
 mannan-MUC1 -- Psiron  
 MAP 30  
 Marinovir -- Phytera  
 MARstem -- Maret  
 MB-015 -- Mochida  
 MBP -- ImmuLogic  
 MCI-028 -- Mitsubishi-Tokyo  
 MCIF -- Human Genome Sciences  
 MDC -- Advanced BioScience -- Akzo  
 Nobel, ICOS  
 MDX 11 -- Medarex  
 MDX 210 -- Medarex  
 MDX 22 -- Medarex  
 MDX 22  
 MDX 240 -- Medarex  
 MDX 33  
 MDX 44 -- Medarex  
 MDX 447 -- Medarex  
 MDX H210 -- Medarex  
 MDX RA -- Houston BioTech., Medarex  
 ME-104 -- Pharmexa  
 Measles vaccine  
 Mecasermin -- Cephalon/Chiron, Chiron  
 MEDI 488 -- MedImmune  
 MEDI 500  
 MEDI 507 -- BioTransplant  
 melanin concentrating hormone --  
 Neurocrine Biosciences  
 melanocortins -- OMRF  
 Melanoma monoclonal antibodies -- Viragen  
 melanoma vaccine -- GlaxoSmithKline,  
 Akzo Nobel, Avant, Aventis Pasteur,  
 Bavarian Nordic, Biovector, CancerVax,  
 Genzyme Molecular Oncology, Humbolt,  
 ImClone Systems, Memorial, NYU, Oxxon  
 Melanoma vaccine Magevac -- Therion  
 memory enhancers -- Scios  
 meningococcal B vaccine -- Chiron  
 meningococcal vaccine -- CAMR  
 Meningococcal vaccine group B conjugate -  
 - North American Vaccine  
 Meningococcal vaccine group B  
 recombinant -- BioChem Vaccines,  
 Microscience  
 Meningococcal vaccine group Y conjugate -  
 - North American Vaccine  
 Meningococcal vaccine groups A B and C  
 conjugate -- North American Vaccine  
 Mepolizumab -- GlaxoSmithKline  
 Metastatin -- EntreMed, Takeda  
 Met-CkB7 -- Human Genome Sciences  
 met-enkephalin -- TNI  
 METH-1 -- Human Genome Sciences  
 methioninase -- AntiCancer  
 Methionine lyase gene therapy --  
 AntiCancer  
 Met-RANTES -- Genexa Biomedical,  
 Seroxo  
 Metreleptin  
 Microtubule inhibitor MAb  
 Immunogen/Abgenix  
 MGDF -- Kirin  
 MGX -- Progenics  
 micrin -- Endocrine  
 microplasmin -- ThromboGenics  
 MIF -- Genetics Institute

FIG. 28R

## 49/498

migration inhibitory factor -- NIH	MAb 45-2D9- -- haematoporphyrin conjugate
Mim CD4.1 -- Xycte Therapies	MAb 4B4
mirostipen -- Human Genome Sciences	MAb 4E3-CPA conjugate -- BCM Oncologia
Mitumomab (BEC-2) -- ImClone Systems, Merck KGaA	MAb 4E3-daunorubicin conjugate
MK 852 -- Merck	MAb 50-6
MLN 1202 (Anti-CCR2 monoclonal antibody) -- Millenium Pharmaceuticals	MAb 50-61A -- Institut Pasteur
Mobenakin -- NIS	MAb 5A8 -- Biogen
molgramostim -- Genetics Institute, Novartis	MAb 791T/36-methotrexate conjugate
monoclonal antibodies -- Abgenix/Celltech, Immusol/ Medarex, Viragen/ Roslin Institute, Cambridge Antibody Tech./Elan	MAb 7c11.e8
MAb 108 --	MAb 7E11 C5-selenocystamine conjugate
MAb 10D5 --	MAb 93KA9 -- Novartis
MAb 14.18-interleukin-2 immunocytokine -- Lexigen	MAb A5B7-cisplatin conjugate -- Biodynamics Research, Pharmacia
MAb 14G2a --	MAb A5B7-I-131
MAb 15A10 --	MAb A7
MAb 170 -- Biomira	MAb A717 -- Exocell
MAb 177Lu CC49 --	MAb A7-zinostatin conjugate
MAb 17F9	MAb ABX-RB2 -- Abgenix
MAb 1D7	MAb ACA 11
MAb 1F7 -- Immune Network	MAb AFP-I-131 -- Immunomedics
MAb 1H10-doxorubicin conjugate	MAb AP1
MAb 26-2F	MAb AZ1
MAb 2A11	MAb B3-LysPE40 conjugate
MAb 2E1 -- RW Johnson	MAb B4 -- United Biomedical
MAb 2F5	MAb B43 Genistein-conjugate
MAb 31.1 -- International BioImmune Systems	MAb B43.13-Tc-99m -- Biomira
MAb 32 -- Cambridge Antibody Tech., Peptech	MAb B43-PAP conjugate
MAb 323A3 -- Centocor	MAb B4G7-gelonin conjugate
MAb 3C5	MAb BCM 43-daunorubicin conjugate -- BCM Oncologia
MAb 3F12	MAb BIS-1
MAb 3F8	MAb BMS 181170 -- Bristol-Myers Squibb
MAb 42/6	MAb BR55-2
MAb 425 -- Merck KGaA	MAb BW494
MAb 447-52D -- Merck Sharp & Dohme	MAb C 242-DM1 conjugate -- ImmunoGen
	MAb C242-PE conjugate
	MAb c30-6
	MAb CA208-cytorhodin-S conjugate -- Hoechst Japan
	MAb CC49 -- Enzo

## 50/498

MAb ch14.18 --	MAb LL2-I-131 -- Immunomedics
MAb CH14.18-GM-CSF fusion protein --	MAb LL2-Y-90
Lexigen	MAb LS2D617 -- Hybritech
MAb chCE7	MAb LYM-1-gelonin conjugate
MAb CI-137 -- AMRAD	MAb LYM-1-I-131
MAb cisplatin conjugate	MAb LYM-1-Y-90
MAb CLB-CD19	MAb LYM-2 -- Peregrine
MAb CLB-CD19v	MAb M195
MAb CLL-1 -- Peregrine	MAb M195-bismuth 213 conjugate --
MAb CLL-1-GM-CSF conjugate	Protein Design Labs
MAb CLL-1-IL-2 conjugate -- Peregrine	MAb M195-gelonin conjugate
MAb CLN IgG -- doxorubicin conjugates	MAB M195-I-131
MAB conjugates -- Tanox	MAB M195-Y-90
MAB D612	MAB MA 33H1 -- Sanofi
MAB Dal B02	MAB MAD11
MAB DC101 -- ImClone	MAB MGB2
MAB EA 1 --	MAB MINT5
MAB EC708 -- Biovation	MAB MK2-23
MAB EP-5C7 -- Protein Design Labs	MAB MOC31 ETA(252-613) conjugate
MAB ERIC-1 -- ICRT	MAB MOC-31-In-111
MAB F105 gene therapy	MAB MOC-31-PE conjugate
MAB FC 2.15	MAB MR6 --
MAB G250 -- Centocor	MAB MRK-16 -- Aventis Pasteur
MAB GA6	MAB MS11G6
MAB GA733	MAB MX-DTPA BrE-3
MAB Gliomab-H -- Viventia Biotech	MAB MY9
MAB HB2-saporin conjugate	MAB Nd2 -- Tosoh
MAB HD 37 --	MAB NG-1 -- Hygeia
MAB HD37-ricin chain-A conjugate	MAB NM01 -- Nissin Food
MAB HNK20 -- Acambis	MAB OC 125
MAB huN901-DM1 conjugate --	MAB OC 125-CMA conjugate
ImmunoGen	MAB OKI-1 -- Ortho-McNeil
MAB I-131 CC49 -- Corixa	MAB OX52 -- Bioproducts for Science
MAB ICO25	MAB PMA5
MAB ICR12-CPG2 conjugate	MAB PR1
MAB ICR-62	MAB prost 30
MAB IRac-ricin A conjugate	MAB R-24
MAB K1	MAB R-24 $\alpha$ Human GD3 -- Celltech
MAB KS1-4-methotrexate conjugate	MAB RFB4-ricin chain A conjugate
MAB L6 -- Bristol-Myers Squibb, Oncogen	MAB RFT5-ricin chain A conjugate
MAB LICO 16-88	MAB SC 1

FIG. 28T

## 51/498

MAb SM-3 -- ICRT  
 MAb SMART 1D10 -- Protein Design Labs  
 MAb SMART ABL 364 -- Novartis  
 MAb SN6f  
 MAb SN6f-deglycosylated ricin A chain conjugate --  
 MAb SN6j  
 MAb SN7-ricin chain A conjugate  
 MAb T101-Y-90 conjugate -- Hybritech  
 MAb T-88 -- Chiron  
 MAb TB94 -- Cancer Immunobiology  
 MAb TEC 11  
 MAb TES-23 -- Chugai  
 MAb TM31 -- Avant  
 MAb TNT-1 -- Cambridge Antibody Tech., Peregrine  
 MAb TNT-3  
 MAb TNT-3 -- IL2 fusion protein --  
 MAb TP3-At-211  
 MAb TP3-PAP conjugate --  
 MAb UJ13A -- ICRT  
 MAb UN3  
 MAb ZME-018-gelonin conjugate  
 MAb-BC2 -- GlaxoSmithKline  
 MAb-DM1 conjugate -- ImmunoGen  
 MAb-ricin-chain-A conjugate -- XOMA  
 MAb-temoporfin conjugates  
 Monopharm C -- Viventia Biotech  
 montepelase -- Eisai  
 montirelin hydrate -- Gruenenthal  
 morotocog alfa -- Genetics Institute  
 Morotocog-alfa -- Pharmacia  
 MP 4  
 MP-121 -- Biopharm  
 MP-52 -- Biopharm  
 MRA -- Chugai  
 MS 28168 -- Mitsui Chemicals, Nihon Schering  
 MSH fusion toxin -- Ligand  
 MSI-99 -- Genæra  
 MT 201 -- Micromet  
 Muc-1 vaccine -- Corixa  
 mucosal tolerance -- Aberdeen  
 mullerian inhibiting subst  
 muplestim -- Genetics Institute, Novartis,  
 DSM Anti-Infectives  
 murine MAb -- KS Biomedix  
 Mutant somatropin -- JCR Pharmaceutical  
 MV 833 -- Toagosei  
 Mycoplasma pulmonis vaccine  
 Mycoprex -- XOMA  
 myeloperoxidase -- Henogen  
 myostatin -- Genetics Institute  
 Nacolomab tafanatox -- Pharmacia  
 Nagrecor -- Scios  
 nagrestipen -- British Biotech  
 NAP-5 -- Corvas Intl.  
 NAPc2 -- Corvas Intl.  
 nartograstim -- Kyowa  
 Natalizumab -- Protein Design Labs  
 Nateplase -- NIH, Nihon Schering  
 nateplase -- Schering AG  
 NBI-3001 -- Neurocrine Biosci.  
 NBI-5788 -- Neurocrine Biosci.  
 NBI-6024 -- Neurocrine Biosci.  
 Nef inhibitors -- BRI  
 Neisseria gonorrhoea vaccine -- Antex Biologics  
 Neomycin B-arginine conjugate  
 Nerelimomab -- Chiron  
 Nerve growth factor -- Amgen -- Chiron, Genentech  
 Nerve growth factor gene therapy  
 nesiritide citrate -- Scios  
 neuregulin-2 -- CeNeS  
 neurocan -- NYU  
 neuronal delivery system -- CAMR  
 Neurophil inhibitory Factor -- Corvas  
 Neuroprotective vaccine -- University of Auckland  
 neurotrophic chimaeras -- Regeneron  
 neurotrophic factor -- NsGene, CereMedix

FIG. 28U

## 52/498

NeuroVax -- Immune Response  
 neurturnin -- Genentech  
 neutral endopeptidase -- Genentech  
 NGF enhancers -- NeuroSearch  
 NHL vaccine -- Large Scale Biology  
 NIP45 -- Boston Life Sciences  
 NKI-B20  
 NM 01 -- Nissin Food  
 NMI-139 -- NitroMed  
 NMMP -- Genetics Institute  
 NN-2211 -- Novo Nordisk  
 Noggin -- Regeneron  
 Nonacog alfa  
 Norelin -- Biostar  
 Norwalk virus vaccine  
 NRLU 10 -- NeoRx  
 NRLU 10 PE -- NeoRx  
 NT-3 -- Regeneron  
 NT-4/5 -- Genentech  
 NU 3056  
 NU 3076  
 NX 1838 -- Gilead Sciences  
 NY ESO-1/CAG-3 antigen -- NIH  
 NYVAC-7 -- Aventis Pasteur  
 NZ-1002 -- Novazyme  
 obesity therapy -- Nobex  
 OC 10426 -- Ontogen  
 OC 144093 -- Ontogen  
 OCIF -- Sankyo  
 Oct-43 -- Otsuka  
 Odulimomab -- Immunotech  
 OK PSA - liposomal  
 OKT3-gamma-1-ala-ala  
 OM 991  
 OM 992  
 Omalizumab -- Genentech  
 oncoimmunin-L -- NIH  
 Oncolysin B -- ImmunoGen  
 Oncolysin CD6 -- ImmunoGen  
 Oncolysin M -- ImmunoGen  
 Oncolysin S -- ImmunoGen  
 Oncophage -- Antigenics  
 Oncostatin M -- Bristol-Myers Squibb  
 OncoVax-CL -- Jenner Biotherapies  
 OncoVax-P -- Jenner Biotherapies  
 onercept -- Yeda  
 onychomycosis vaccine -- Boehringer  
 Ingelheim  
 opebecan -- XOMA  
 opioids -- Arizona  
 Oprelvekin -- Genetics Institute  
 Oregovomab -- AltaRex  
 Org-33408 b-- Akzo Nobel  
 Orlip DP -- EpiCept  
 oryzacystatin  
 OSA peptides -- GenSci Regeneration  
 osteoblast-cadherin GF -- Pharis  
 Osteocalcin-thymidine kinase gene therapy  
 osteogenic protein -- Curis  
 osteopontin -- OraPharma  
 osteoporosis peptides -- Integra, Telios  
 osteoprotegerin -- Amgen, SnowBrand  
 otitis media vaccines -- Antex Biologics  
 ovarian cancer -- University of Alabama  
 OX40-IgG fusion protein -- Cantab, Xenova  
 P 246 -- Diatide  
 P 30 -- Alfacell  
 p1025 -- Active Biotech  
 P-113<sup>A</sup> -- Demegen  
 P-16 peptide -- Transition Therapeutics  
 p43 -- Ramot  
 P-50 peptide -- Transition Therapeutics  
 p53 + RAS vaccine -- NIH, NCI  
 PACAP(1-27) analogue  
 paediatric vaccines -- Chiron  
 Pafase -- ICOS  
 PAGE-4 plasmid DNA -- IDEC  
 PAI-2 -- Biotech Australia, Human  
 Therapeutics  
 Palifermin (keratinocyte growth factor) --  
 Amgen  
 Palivizumab -- MedImmune

FIG. 28V

53/498

PAM 4 -- Merck  
 pamiteplase -- Yamanouchi  
 pancreatin, Minitabs -- Eurand  
 Pangen -- Fournier  
 Pantarin -- Selective Genetics  
 Parainfluenza virus vaccine -- Pharmacia,  
 Pierre Fabre  
 paraoxanase -- Esperion  
 parathyroid hormone -- Abiogen, Korea  
 Green Cross  
 Parathyroid hormone (1-34) --  
 Chugai/Suntory  
 Parkinson's disease gene therapy -- Cell  
 Genesys/ Ceregene  
 Parvovirus vaccine -- MedImmune  
 PCP-Scan -- Immunomedics  
 PDGF -- Chiron  
 PDGF cocktail -- Theratechnologies  
 peanut allergy therapy -- Dynavax  
 PEG anti-ICAM MAb -- Boehringer  
 Ingelheim  
 PEG asparaginase -- Enzon  
 PEG glucocerebrosidase  
 PEG hirudin -- Knoll  
 PEG interferon-alpha-2a -- Roche  
 PEG interferon-alpha-2b + ribavirin --  
 Biogen, Enzon, ICN Pharmaceuticals,  
 Schering-Plough  
 PEG MAb A5B7 --  
 Pegacaristim -- Amgen -- Kirin Brewery --  
 ZymoGenetics  
 Pegaldesleukin -- Research Corp  
 pegaspargase -- Enzon  
 pegfilgrastim -- Amgen  
 PEG-interferon Alpha -- Viragen  
 PEG-interferon Alpha 2A -- Hoffman La-  
 Roche  
 PEG-interferon Alpha 2B -- Schering-  
 Plough  
 PEG-r-hirudin -- Abbott  
 PEG-rHuMGDF -- Amgen  
 PEG-uricase -- Mountain View  
 Pegvisomant -- Genentech  
 PEGylated proteins, PolyMASC -- Valentis  
 PEGylated recombinant native human leptin  
 -- Roche  
 Pentumomab  
 Penetratin -- Cyclacel  
 Pepscan -- Antisoma  
 peptide G -- Peptech, ICRT  
 peptide vaccine -- NIH ,NCI  
 Pexelizumab  
 pexiganan acetate -- Genaera  
 Pharmaprojects No. 3179 -- NYU  
 Pharmaprojects No. 3390 -- Ernest Orlando  
 Pharmaprojects No. 3417 -- Sumitomo  
 Pharmaprojects No. 3777 -- Acambis  
 Pharmaprojects No. 4209 -- XOMA  
 Pharmaprojects No. 4349 -- Baxter Intl.  
 Pharmaprojects No. 4651  
 Pharmaprojects No. 4915 -- Avanir  
 Pharmaprojects No. 5156 -- Rhizogenics  
 Pharmaprojects No. 5200 -- Pfizer  
 Pharmaprojects No. 5215 -- Origene  
 Pharmaprojects No. 5216 -- Origene  
 Pharmaprojects No. 5218 -- Origene  
 Pharmaprojects No. 5267 -- ML  
 Laboratories  
 Pharmaprojects No. 5373 -- MorphoSys  
 Pharmaprojects No. 5493 -- Metabolex  
 Pharmaprojects No. 5707 -- Genentech  
 Pharmaprojects No. 5728 -- Autogen  
 Pharmaprojects No. 5733 -- BioMarin  
 Pharmaprojects No. 5757 -- NIH  
 Pharmaprojects No. 5765 -- Gryphon  
 Pharmaprojects No. 5830 -- AntiCancer  
 Pharmaprojects No. 5839 -- Dyax  
 Pharmaprojects No. 5849 -- Johnson &  
 Johnson  
 Pharmaprojects No. 5860 -- Mitsubishi-  
 Tokyo

FIG. 28W

## 54/498

- Pharmaprojects No. 5869 -- Oxford GlycoSciences
- Pharmaprojects No. 5883 -- Asahi Brewery
- Pharmaprojects No. 5947 -- StressGen
- Pharmaprojects No. 5961 -- Theratechnologies
- Pharmaprojects No. 5962 -- NIH
- Pharmaprojects No. 5966 -- NIH
- Pharmaprojects No. 5994 -- Pharming
- Pharmaprojects No. 5995 -- Pharming
- Pharmaprojects No. 6023 -- IMMUCON
- Pharmaprojects No. 6063 -- Cytoclonal
- Pharmaprojects No. 6073 -- SIDDICO
- Pharmaprojects No. 6115 -- Genzyme
- Pharmaprojects No. 6227 -- NIH
- Pharmaprojects No. 6230 -- NIH
- Pharmaprojects No. 6236 -- NIH
- Pharmaprojects No. 6243 -- NIH
- Pharmaprojects No. 6244 -- NIH
- Pharmaprojects No. 6281 -- Senetek
- Pharmaprojects No. 6365 -- NIH
- Pharmaprojects No. 6368 -- NIH
- Pharmaprojects No. 6373 -- NIH
- Pharmaprojects No. 6408 -- Pan Pacific
- Pharmaprojects No. 6410 -- Athersys
- Pharmaprojects No. 6421 -- Oxford GlycoSciences
- Pharmaprojects No. 6522 -- Maxygen
- Pharmaprojects No. 6523 -- Pharis
- Pharmaprojects No. 6538 -- Maxygen
- Pharmaprojects No. 6554 -- APALEXO
- Pharmaprojects No. 6560 -- Ardana
- Pharmaprojects No. 6562 -- Bayer
- Pharmaprojects No. 6569 -- Eos
- Phenoxazine
- Phenylase -- Ibbex
- Pigment epithelium derived factor -- plasminogen activator inhibitor-1, recombinant -- DuPont Pharmaceuticals
- Plasminogen activators -- Abbott Laboratories, American Home Products, Boehringer Mannheim, Chiron Corporation, DuPont Pharmaceuticals, Eli Lilly, Shionogi, Genentech, Genetics Institute, GlaxoSmithKline, Hemispherx Biopharma, Merck & Co, Novartis, Pharmacia Corporation, Wakamoto, Yeda
- plasminogen-related peptides -- Bio-Tech. General/MGH
- platelet factor 4 -- RepliGen
- Platelet-derived growth factor -- Amgen -- ZymoGenetics
- plusonemin -- Hayashibara
- PMD-2850 -- Protherics
- Pneumococcal vaccine -- Antex Biologics, Aventis Pasteur
- Pneumococcal vaccine intranasal -- BioChem Vaccines/Biovector
- PR1A3
- PR-39
- pralmorelin -- Kaken
- Pretarget-Lymphoma -- NeoRx
- Priliximab -- Centocor
- PRO 140 -- Progenics
- PRO 2000 -- Procept
- PRO 367 -- Progenics
- PRO 542 -- Progenics
- pro-Apo A-I -- Esperion
- prolactin -- Genzyme
- Prosaptide TX14(A) -- Bio-Tech. General
- prostate cancer antibodies -- Immunex, UroCor
- prostate cancer antibody therapy -- Genentech/UroGenesys, Genotherapeutics
- prostate cancer immunotherapeutics -- The PSMA Development Company
- prostate cancer vaccine -- Aventis Pasteur, Zonagen, Corixa, Dendreon, Jenner Biotherapies, Therion Biologics

FIG. 28X

## 55/498

prostate-specific antigen -- Entremed  
 protein A -- RepliGen  
 protein adhesives -- Enzon  
 protein C -- Baxter Intl., PPL Therapeutics,  
     ZymoGenetics  
 protein C activator -- Gilead Sciences  
 protein kinase R antagonists -- NIH  
 protirelin -- Takeda  
 protocadherin 2 -- Caprion  
 Pro-urokinase -- Abbott, Bristol-Myers  
     Squibb, Dainippon, Tosoh -- Welfide  
 P-selectin glycoprotein ligand-1 -- Genetics  
     Institute  
 pseudomonal infections -- InterMune  
 Pseudomonas vaccine -- CytoVax  
 PSGL-Ig -- American Home Products  
 PSP-94 -- Procyon  
 PTH 1-34 -- Nobex  
 Quilimmune-M -- Antigenics  
 R 744 -- Roche  
 R 101933  
 R 125224 -- Sankyo  
 RA therapy -- Cardion  
 Rabies vaccine recombinant -- Aventis  
     Pasteur, BioChem Vaccines, Kaketsuken  
     Pharmaceuticals  
 RadioTheraCIM -- YM BioSciences  
 Ramot project No. 1315 -- Ramot  
 Ramot project No. K-734A -- Ramot  
 Ramot project No. K-734B -- Ramot  
 Ranibizumab (Anti-VEGF fragment) --  
     Genentech  
 RANK -- Immunex  
 ranpirinase -- Alfacell  
 ranpirinase-anti-CD22 MAb -- Alfacell  
 RANTES inhibitor -- Milan  
 RAPID drug delivery systems -- ARIAD  
 rasburicase -- Sanofi  
 rBPI-21, topical -- XOMA  
 RC 529 -- Corixa  
 rCFTR -- Genzyme Transgenics

RD 62198  
 rDnase -- Genentech  
 RDP-58 -- SangStat  
 ReceptTox-Fce -- Keryx  
 ReceptTox-GnRH -- Keryx, MTR  
     Technologies  
 ReceptTox-MBP -- Keryx, MTR  
     Technologies  
 recFSH -- Akzo Nobel, Organon  
 REGA 3G12  
 Regavirumab -- Teijin  
 relaxin -- Connetics Corp  
 Renal cancer vaccine -- MacroPharm  
 repifermin -- Human Genome Sciences  
 Respiratory syncytial virus PFP-2 vaccine --  
     Wyeth-Lederle  
 Respiratory syncytial virus vaccine --  
     GlaxoSmithKline, Pharmacia, Pierre Fabre  
 Respiratory syncytial virus vaccine  
     inactivated  
 Respiratory syncytial virus-parainfluenza  
     virus vaccine -- Aventis Pasteur,  
     Pharmacia  
 Reteplase -- Boehringer Mannheim,  
     Hoffman La-Roche  
 Retropep -- Retroscreen  
 RFB4 (dsFv) PE38  
 RFI 641 -- American Home Products  
 RFTS -- UAB Research Foundation  
 RG 12986 -- Aventis Pasteur  
 RG 83852 -- Aventis Pasteur  
 RG-1059 -- RepliGen  
 rGCR -- NIH  
 rGLP-1 -- Restoragen  
 rGRF -- Restoragen  
 rh Insulin -- Eli Lilly  
 RHAMM targeting peptides -- Cangene  
 rHb1.1 -- Baxter Intl.  
 rhCC10 -- Claragen  
 rhCG -- SeroNo  
 Rheumatoid arthritis gene therapy

FIG. 28Y



## 56/498

Rheumatoid arthritis vaccine -- Veterans

Affairs Medical Center

rhLH -- SeroNo

Ribozyme gene therapy -- Genset

Rickettsial vaccine recombinant

RIGScan CR -- Neoprobe

RIP-3 -- Rigel

Rituximab -- Genentech

RK-0202 -- RxKinetix

RLT peptide -- Esperion

rM/NEI -- IVAX

rmCRP -- Immtech

RN-1001 -- Renovo

RN-3 -- Renovo

RNase conjugate -- Immunomedics

RO 631908 -- Roche

Rotavirus vaccine -- Merck

RP 431 -- DuPont Pharmaceuticals

RP-128 -- Resolution

RPE65 gene therapy --

RPR 110173 -- Aventis Pasteur

RPR 115135 -- Aventis Pasteur

RPR 116258A -- Aventis Pasteur

rPSGL-Ig -- American Home Products

r-SPC surfactant -- Byk Gulden

RSV antibody -- Medimmune

Ruplizumab -- Biogen

rV-HER-2/neu -- Therion Biologics

SA 1042 -- Sankyo

sacrosidase -- Orphan Medical

Sant 7

Sargramostim -- Immunex

saruplase -- Gruenenthal

Satumomab -- Cytogen

SB 1 -- COR Therapeutics

SB 207448 -- GlaxoSmithKline

SB 208651 -- GlaxoSmithKline

SB 240683 -- GlaxoSmithKline

SB 249415 -- GlaxoSmithKline

SB 249417 -- GlaxoSmithKline

SB 6 -- COR Therapeutics

SB RA 31012 --

SC 56929 -- Pharmacia

SCA binding proteins -- Curis, Enzon

scFv(14E1)-ETA Berlex Laboratories,

Schering AG

ScFv(FRP5)-ETA --

ScFv6C6-PE40 --

SCH 55700 -- Celltech

Schistosomiasis vaccine -- Glaxo

Wellcome/Medeva, Brazil

SCPF -- Advanced Tissue Sciences

scuPA-suPAR complex -- Hadasit

SD-9427 -- Pharmacia

SDF-1 -- Ono

SDZ 215918 -- Novartis

SDZ 280125 -- Novartis

SDZ 89104 -- Novartis

SDZ ABL 364 -- Novartis

SDZ MMA 383 -- Novartis

Secretin -- Ferring, Repligen

serine protease inhbs -- Pharis

sermorelin acetate -- SeroNo

SERP-1 -- Viron

sertenef -- Dainippon

serum albumin, Recombinant human --

Aventis Behring

serum-derived factor -- Hadasit

Sevirumab -- Novartis

SGN 14 -- Seattle Genetics

SGN 15 -- Seattle Genetics

SGN 17/19 -- Seattle Genetics

SGN 30 -- Seattle Genetics

SGN-10 -- Seattle Genetics

SGN-11 -- Seattle Genetics

SH 306 -- DuPont Pharmaceuticals

Shanvac-B -- Shantha

Shigella flexneri vaccine -- Avant, Acambis,

Novavax

Shigella sonnei vaccine --

sICAM-1 -- Boehringer Ingelheim

Silteplase -- Genzyme

**FIG. 28Z**

## 57/498

SIV vaccine -- Endocon, Institut Pasteur  
 SK 896 -- Sanwa Kagaku Kenkyusho  
 SK-827 -- Sanwa Kagaku Kenkyusho  
 Skeletex -- CellFactors  
 SKF 106160 -- GlaxoSmithKline  
 S-nitroso-AR545C --  
 SNTP -- Active Biotech  
 somatomedin-1 -- GroPep, Mitsubishi-  
 Tokyo, NIH  
 somatomedin-1 carrier protein -- Insmed  
 somatostatin -- Ferring  
 Somatotropin/  
 Human Growth Hormone -- Bio-Tech.  
 General, Eli Lilly  
 somatropin -- Bio-Tech. General, Alkermes,  
 ProLease, Aventis Behring, Biovector,  
 Cangene, Dong-A, Eli Lilly, Emisphere,  
 Enact, Genentech, Genzyme Transgenics,  
 Grandis/InfiMed, CSL, InfiMed, MacroMed,  
 Novartis, Novo Nordisk, Pharmacia  
 Serono, TranXenoGen  
 somatropin derivative -- Schering AG  
 somatropin, AIR -- Eli Lilly  
 Somatropin, inhaled -- Eli Lilly/Alkermes  
 somatropin, Kabi -- Pharmacia  
 somatropin, Orasome -- Novo Nordisk  
 Sonermin -- Dainippon Pharmaceutical  
 SP(V5.2)C -- Supertek  
 SPf66  
 sphingomyelinase -- Genzyme  
 SR 29001 -- Sanofi  
 SR 41476 -- Sanofi  
 SR-29001 -- Sanofi  
 SS1(dsFV)-PE38 -- NeoPharm  
 $\beta$ 2 microglobulin -- Avidex  
 $\beta$ 2-microglobulin fusion proteins -- NIH  
 $\beta$ -amyloid peptides -- CeNeS  
 $\beta$ -defensin -- Pharis  
 Staphylococcus aureus infections --  
 Inhibibex/ZLB  
 Staphylococcus aureus vaccine conjugate --  
 Nabi  
 Staphylococcus therapy -- Tripep  
 Staphylokinase -- Biovation, Prothera,  
 Thrombogenetics  
 Streptococcal A vaccine -- M6  
 Pharmaceuticals, North American Vaccine  
 Streptococcal B vaccine -- Microscience  
 Streptococcal B vaccine recombinant --  
 Biochem Vaccines  
 Streptococcus pyogenes vaccine  
 STRL-33 -- NIH  
 Subalin -- SRC VB VECTOR  
 SUIIS -- United Biomedical  
 SUIIS-LHRH -- United Biomedical  
 SUN-E3001 -- Suntory  
 super high affinity monoclonal antibodies --  
 YM BioSciences  
 Superoxide dismutase -- Chiron, Enzon,  
 Ube Industries, Bio-Tech, Yeda  
 superoxide dismutase-2 -- OXIS  
 suppressin -- UAB Research Foundation  
 SY-161-P5 -- ThromboGenics  
 SY-162 -- ThromboGenics  
 Systemic lupus erythematosus vaccine --  
 MedClone/VivoRx  
 T cell receptor peptides -- Xoma  
 T cell receptor peptide vaccine  
 T4N5 liposomes -- AGI Dermatics  
 TACI, soluble -- ZymoGenetics  
 targeted apoptosis -- Antisoma  
 tasonermin -- Boehringer Ingelheim  
 TASP  
 TASP-V  
 Tat peptide analogues -- NIH  
 TBP I -- Yeda  
 TBP II  
 TBV25H -- NIH  
 Tc 99m ior cea1 -- Center of Molecular  
 Immunology  
 Tc 99m P 748 -- Diatide

FIG. 28AA

## 58/498

Tc 99m votumab -- Intracell  
 Tc-99m rh-Annexin V -- Theseus Imaging  
 teceleukin -- Biogen  
 tenecteplase -- Genentech  
 Teriparatide -- Armour Pharmaceuticals,  
     Asahi Kasei, Eli Lilly  
 terlipressin -- Ferring  
 testisin -- AMRAD  
 Tetra fibrin -- Roche  
 TFPI -- EntrelMed  
 tgD-IL-2 -- Takeda  
 TGF-Alpha -- ZymoGenetics  
 TGF- $\beta$  -- Kolon  
 TGF- $\beta$ 2 -- Insmed  
 TGF- $\beta$ 3 -- OSI  
 Thalassemia gene therapy -- Crucell  
 TheraCIM-h-R3 -- Center of Molecular  
     Immunology, YM BioSciences  
 Theradigm-HBV -- Epimmune  
 Theradigm-HPV -- Epimmune  
 Theradigm-malaria -- Epimmune  
 Theradigm-melanoma -- Epimmune  
 TheraFab -- Antisoma  
 ThGRF 1-29 -- Theratechnologies  
 ThGRF 1-44 -- Theratechnologies  
 Thrombin receptor activating peptide --  
     Abbott  
 thrombomodulin -- Iowa, Novocastra  
 Thrombopoietin -- Dragon Pharmaceuticals,  
     Genentech  
 thrombopoietin, Pliva -- Recepton  
 Thrombospondin 2 --  
 thrombostatin -- Thromgen  
 thymalfasin -- SciClone  
 thymocartin -- Gedeon Richter  
 thymosin Alpha1 -- NIH  
 thyroid stimulating hormone -- Genzyme  
 tICAM-1 -- Bayer  
 Tick anticoagulant peptide -- Merck  
 TIF -- Xoma  
 Tifacogin -- Chiron, NIS, Pharmacia  
 Tissue factor -- Genentech  
 Tissue factor pathway inhibitor  
 TJN-135 -- Tsumura  
 TM 27 -- Avant  
 TM 29 -- Avant  
 TMC-151 -- Tanabe Seiyaku  
 TNF tumour necrosis factor -- Asahi Kasei  
 TNF Alpha -- CytImmune  
 TNF antibody -- Johnson & Johnson  
 TNF binding protein -- Amgen  
 TNF degradation product -- Oncotech  
 TNF receptor -- Immunex  
 TNF receptor 1, soluble -- Amgen  
 TNF Tumour necrosis factor-alpha -- Asahi  
     Kasei, Genentech, Mochida  
 TNF-Alpha inhibitor -- Tripep  
 TNFR:Fc gene therapy -- Targeted Genetics  
 TNF-SAM2  
 Tolerimab -- Innogenetics  
 Toxoplasma gondii vaccine --  
     GlaxoSmithKline  
 TP 9201 -- Telios  
 TP10 -- Avant  
 TP20 -- Avant  
 tPA -- Centocor  
 trafermin -- Scios  
 TRAIL/Apo2L -- Immunex  
 TRAIL-R1 MAb -- Cambridge Antibody  
     Technologies  
 transferrin-binding proteins -- CAMR  
 Transforming growth factor-beta-1 --  
     Genentech  
 transport protein -- Genesis  
 Trastuzumab -- Genentech  
 TRH -- Ferring  
 Triabin -- Schering AG  
 Triconal  
 Triflavin  
 troponin I -- Boston Life Sciences  
 TRP-2<sup>A</sup> -- NIH  
 trypsin inhibitor -- Mochida

FIG. 28BB

## 59/498

<p>TSP-1 gene therapy –  TT-232  TTS-CD2 -- Active Biotech  Tuberculosis vaccine -- Aventis Pasteur,  Genesis  Tumor Targeted Superantigens -- Active  Biotech -- Pharmacia  tumour vaccines -- PhotoCure  tumour-activated prodrug antibody  conjugates -- Millennium/ImmunoGen  tumstatin -- ILEX  Tuvirumab -- Novartis  TV-4710 – Teva  TWEAK receptor -- Immunex  TXU-PAP  TY-10721 – TOA Eiyo  Type I diabetes vaccine -- Research Corp  Typhoid vaccine CVD 908  U 143677 -- Pharmacia  U 81749 -- Pharmacia  UA 1248 -- Arizona  UGIF -- Sheffield  UIC 2  UK 101  UK-279276 -- Corvas Intl.  urodilatin -- Pharis  urofollitrophin -- Serono  Urokinase -- Abbott  uteroferin-- Pepgen  V 20 -- GLYCODesign  V2 vasopressin receptor gene therapy  vaccines -- Active Biotech  Varicella zoster glycoprotein vaccine --  Research Corporation Technologies  Varicella zoster virus vaccine live -- Cantab  Pharmaceuticals  Vascular endothelial growth factor --  Genentech, University of California</p>	<p>Vascular endothelial growth factors -- R&amp;D  Systems  vascular targeting agents -- Peregrine  vasopermeation enhancement agents --  Peregrine  vasostatin -- NIH  VCL -- Bio-Tech. General  VEGF -- Genentech, Scios  VEGF inhibitor -- Chugai  VEGF-2 -- Human Genome Sciences  VEGF-Trap -- Regeneron  viscumin, recombinant -- Madaus  Vitaxin  Vitrax -- ISTA Pharmaceuticals  West Nile virus vaccine -- Bavarian Nordic  WP 652  WT1 vaccine -- Corixa  WX-293 -- Wilex BioTech.  WX-360 -- Wilex BioTech.  WX-UK1 -- Wilex BioTech.  XMP-500 -- XOMA  XomaZyme-791 -- XOMA  XTL 001 -- XTL Biopharmaceuticals  XTL 002 -- XTL Biopharmaceuticals  yeast delivery system -- GlobelImmune  Yersinia pestis vaccine  YIGSR-Stealth -- Johnson &amp; Johnson  Yissum Project No. D-0460 -- Yissum  YM 207 -- Yamanouchi  YM 337 -- Protein Design Labs  Yttrium-90 labelled biotin  Yttrium-90-labeled anti-CEA MAb T84.66 --  ZD 0490 -- AstraZeneca  ziconotide -- Elan  ZK 157138 -- Berlex Laboratories  Zolimomab aritox  Zorcell -- Immune Response  ZRXL peptides -- Novartis</p>
---	---

FIG. 28CC